

NuAmpsTM User's Guide

Compumedics Neuroscan

7850 Paseo Del Norte
El Paso, TX 79912 USA

Voice :915-845-5600
Fax: 915-845-2965
Internet: sales@neuro.com
techsup@neuro.com
www.neuro.com



For Technical Support.....

If you have any questions or problems, contact Technical Support through any of the following routes.

If you live outside the USA or Canada, and purchased your system through one of our international distributors, please contact the **distributor** first, especially if your system is under warranty. In all other cases, please use **techsup@neuro.com**, or see the other Support options on our web site (<http://www.neuro.com>). Or, if you live in the USA or Canada, please call **1-800 474-7875**. International callers should use **915-845-5600, ext 324**.

For Sales related questions, please contact your local distributor, or contact us at **sales@neuro.com**.

Copyright © 2003 - Compumedics Neuroscan.

All rights reserved. Printed in the United States of America. No part of this manual may be used or reproduced in any form or by any means, or stored in a database or retrieval system, without prior written permission of the company. Making copies of any part of this document for any purpose other than your own personal use is a violation of United States copyright laws. For information, contact Neuroscan.

Compumedics Neuroscan

7850 Paseo Del Norte

El Paso, TX 79912

Telephone: 915-845-5600

Fax: 915-845-2965

E-mail: techsup@neuro.com

Web: www.neuro.com

Table of Contents

Device Classification	5
Introduction	13
Cascaded NuAmps Systems	14
Installing NuAmps	14
1. Unpacking the NuAmps	15
Cascaded NuAmps Systems	17
2. Connecting the Isolation Transformer	18
3. Install the SCAN Software	18
4. Install the NuAmps Hardware	19
Cascaded NuAmps Systems	20
Driver Verification	22
How many channels are in Cascaded NuAmps Systems?	22
5. Configure SCAN for NuAmps	23
Cascaded NuAmps Systems	24
Differences Between Single and Cascaded Systems	25
Acquisition from Cascaded NuAmps	28
Tuning Performance of Windows XP Systems and Cascaded NuAmps	28
Safety and Quality Warnings	30
6. Quick System Test	30
The Amplifiers	31
Software Control Of NuAmps	32
NuAmps Setup Files	33
Enter Number of Channels	34

Enter Channel Labels.....	36
Channel Layout.....	38
NuAmps Setup Screen.....	39
Reference.....	40
Amplifiers.....	41
Trigger Test.....	43
Creating a Bipolar Montage.....	43
DC Electrode Considerations.....	46
Notch Filters and Noise.....	49
Calibration.....	51
Impedance.....	51
Subject Ground.....	53
Electrode Safety.....	53
FAQs.....	53
Troubleshooting.....	56
Appendix A: Trigger Port Interface.....	61
Appendix B: NuAmps Electrode Starter Kit.....	67
Appendix C: PocketTrace Signal Generator.....	78
Appendix D: Cap Connector Pinout.....	91
Appendix E: Impedance Measurement.....	92

Device Classification



ATTENTION: CONSULT ACCOMPANYING DOCUMENTS BEFORE USING

The NuAmps™ Model 7181 amplifiers are a computer-powered instrument designed to meet the applicable requirements of IEC601-1:1988. NuAmps should be used only according to the manufacturer's instructions. Replacement parts and accessories may be obtained from the manufacturer.

Manufacturer: Compumedics Neuroscan
7850 Paseo Del Norte
El Paso, TX, USA 79912
Phone: 915-845-5600
Fax: 915-845-2965
e-mail: sales@neuro.com
techsup@neuro.com
Web: www.neuro.com

This equipment has been tested and found to comply with the limits and requirements for a Class II (portable version) and I (stationary version), Type BF device per EN60601-1-2 and EN55011. These limits and requirements are designed to provide reasonable protection under conditions of normal use from interference with and by other devices. There is, however, no guarantee that interference will not result from operation of this device in proximity or connected to some

other device. If interference occurs, the user or operator is encouraged to try and correct the interference by one or more of the following measures:

- 1) Change the orientation of the two devices relative to one another.
- 2) Increase the separation between the two devices.
- 3) Check the power source and grounding for the two devices.
- 4) Twist amplifier USB cable with photic stimulator serial and power cables
- 5) Consult the dealer, Neuroscan Technical Support, or an experienced technician for help. A slight degradation of performance may be experienced when the NuAmps amplifier is exposed to an RF environment.

The safety and electromagnetic compatibility of this system was tested with the following accessories, parts, and associated devices. The user or operator is cautioned to ensure that when using accessories, parts, or associated devices other than those listed, that the safety and electromagnetic compatibility of the system is maintained.

Portable system:

- 1) Neuroscan NuAmps amplifier, Model 7181
- 2) Neuroscan NuAmps photic stimulator, Model 7097
- 3) IBM® laptop computer, Model ThinkPad® 390E (CE Marked) or
DELL® laptop computer, Model Inspiron™ 3800 (CE Marked) or
other (equivalent and CE Marked) models
- 4) IPSI 240VAC/15-20VDC medical grade power supply, Model PMP55-13 (CE Marked)
- 5) AULT 240VAC/5VDC medical grade power supply, Model SW170 (CE Marked)
- 6) BROTHER® printer, Model MP-21C/CDX (CE Marked) or
other (equivalent and CE Marked) models

Stationary (Lab) system:

- 1) Neuroscan NuAmps amplifier, Model 7181

- 2) Neuroscan NuAmps photic stimulator, Model 7097
- 3) MICRON desktop computer, Model ClientPro® PIII650 (CE Marked) or other (equivalent and CE Marked) models
- 4) NEC LCD display, Model MultiSync 1525V (CE Marked) or other (equivalent and CE Marked) models
- 5) Toroid 240VAC isolation power supply
- 6) HP printer, Model 952C (CE Marked) or other (equivalent and CE Marked) models

Classification per IEC601-1:1988

The device is ordinary equipment not protected against ingress of water and should not be used in the presence of any spilled liquids. It is not designed to be suitable for use in the presence of a flammable anesthetic mixture of air and oxygen or nitrous oxide. The device is capable of continuous operation.

Class and degree of protection against electrical shock is Class II (portable version) or I (stationary version), Type BF.

Technical Description

NuAmps amplifier, Model 7181:

Input:	5VDC (from USB cable)
Fuses:	None
Weight:	~1.4 lbs (695 grams)
Dimensions:	7.8"x5.9"x1.6" (198x151x40 mm)

NuAmps photic stimulator, Model 7097:

Input:	5VDC (from AULT power supply)
Fuses:	None

Weight: ~0.4 lbs (181 grams)
Dimensions: 4.5"x3.1"x1.3" (114.3x78.2x31.8 mm)

Shipping and Storage Maximum Limits

-10° C to + 50° C, 10% to 100% Non-condensing RH, 500 hPa to 1060 hPa. After unpacking, allow devices to adjust to room temperature for at least two hours prior to interconnection and application of power.

Warnings and Precautions

Instructions

Read instructions before operating the device.

Symbols


The following symbols are found on the NuAmps amplifiers:




CAUTION: Read instructions before use!



Subject connection is isolated from power mains.

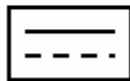

Volts DC

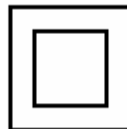

Double Isolation



Manufacturing date

The NuAmps amplifiers contains no user serviceable parts and so are marked “CAUTION: Refer servicing to qualified personnel”.

This device is not equipped with appropriate alarms required for use in monitoring clinical parameters of a patient where it is necessary to alert the user of situations, which could lead to death or severe deterioration of the patient’s state of health.

Manufacturer designs NuAmps amplifier Sync1 and Sync2 ports for connection only to Neuroscan Labs approved accessories. Please contact Neuroscan Labs before connecting a device to these ports.

Environment

The NuAmps amplifiers are designed for usage in a laboratory or clinical environment. Extremes of humidity, temperature, or pressure should be avoided. The device should not be used in a location where contact with liquids is possible, and if liquids are spilled on or in the area of the device, it should not be used until it can be ensured that the fluid or its residue will not affect device operation. Questions should be directed to the manufacturer or its representatives.

Cleaning Instructions

The NuAmps enclosure may be cleaned with a damp sponge or cloth and mild nonabrasive cleanser. Take care to ensure that liquid does not spill in or on the device. Do not use abrasives or detergents.

Sterilization and Cleaning of Patient-Contact Parts

Patient electrodes are supplied as part of the NuAmps system. Follow their manufacturer’s instructions for sterilization and cleaning of these parts. Some devices are designed for onetime use only, and no attempt should be made to reuse them, whether sterilization has been attempted or not. Contact Neuroscan Labs technical support if you have questions about sterilization or cleaning of the NuAmps amplifier or electrodes to be used with the device.

Repair

There are no user serviceable parts neither fuses in the NuAmps amplifier, which is indicated on the back panel label as “CAUTION: Refer servicing to qualified personnel”. Contact your dealer or Neuroscan Labs Technical Support if you believe the NuAmps system is in need of repair.

Maintenance

Neuroscan suggests that the earth and patient leakage currents be tested at least once per year to ensure continued safe use of the device. Also at least once per year visually inspect the device, including power cord. Replace any worn or frayed cables, and contact your dealer or Neuroscan Technical Support if you have concerns about what you see. This inspection interval may be shortened for devices that are moved often or experience unusually heavy use. No calibration of device is required.

Installation Precaution

Proper grounding is important for continued safe use of your NuAmps system. Ensure that the outlet supplying power to your NuAmps system is grounded, and that the power cords supplied with your system are used. Other devices in the same patient area should be at the same ground potential, and should preferably use the same branch circuit. See the Hardware and Software Installation directions below for more details.

Power Source Characteristics

The NuAmps device is designed, produced, and tested to ensure reliable operation when connected to power systems having normal variability. If you believe that your power system may experience excessive noise or variability, Neuroscan recommends use of a power conditioner.

Interconnection with Other Devices

Care should be taken when multiple devices are connected to a patient, or when devices are connected together. Leakage currents for individual devices may sum to values higher than expected for single devices. In particular care should be taken when connecting Information Technology (computer) equipment to Medical equipment. Allowable leakage current levels for IT equipment are higher than for Medical equipment.

Use With HF Surgical Equipment

This device does not contain protection against burning of the patient when used with high frequency (HF) surgical equipment. Neuroscan recommends that the NuAmps device not be connected to the patient during use of HF surgical equipment.

Electrode Safety

The NuAmps amplifier inputs and attached electrodes are Type BF, which means in part that they are not connected to Earth Ground or Chassis Ground. Maintain this separation from Earth Ground by ensuring that the electrodes and any conductive parts of their connectors do not touch conductive parts, including the system enclosure or other grounded devices.

NO PART OF THIS PUBLICATION MAY BE REPRODUCED, STORED IN A RETRIEVAL SYSTEM, OR TRANSMITTED IN ANY FORM OR BY ANY MEANS, MECHANICAL, ELECTRONIC, PHOTOCOPYING, RECORDING, OR OTHERWISE, WITHOUT WRITTEN PERMISSION OF NEUROSCAN.

NO LICENSE, EXPRESS OR IMPLIED, IS GRANTED WITH RESPECT TO ANY OF THE TECHNOLOGY DESCRIBED IN THIS MANUAL. NEUROSCAN RETAINS ALL INTELLECTUAL PROPERTY RIGHTS ASSOCIATED WITH THE TECHNOLOGY DESCRIBED IN THIS MANUAL. THIS MANUAL IS INTENDED TO

ASSIST APPLICATION DEVELOPERS IN THE USE OF APPLICATIONS FOR NUAMPS SYSTEM FROM NEUROSCAN.

EVEN THOUGH NEUROSCAN HAS REVIEWED THIS MANUAL, NEUROSCAN MAKES NO WARRANTY OR REPRESENTATION, EITHER EXPRESS OR IMPLIED, WITH RESPECT TO THIS MANUAL, ITS QUALITY, ACCURACY, MERCHANTABILITY, OR FITNESS FOR A PARTICULAR PURPOSE.

IN NO EVENT WILL NEUROSCAN BE LIABLE FOR DIRECT, INDIRECT, SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES RESULTING FROM ANY DEFECT OR INACCURACY IN THIS MANUAL, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

NuAmps is a trademark of Compumedics Neuroscan.

Windows, Windows 98, Windows 2000, Windows XP, Velcro, IBM, ThinkPad, Dell, Inspiron, Brother, Micro, ClientPro, NEC, MultiSync, HP, Ault, Dove and Ivory are trademarks or registered trademarks of their respective owners.

Introduction

The Neuroscan NuAmps are DC amplifiers designed to record a wide variety of multichannel neurophysiological signals, such as EEG, ECG, EOG, and EMG. The NuAmps are intended for the researcher who needs a small portable amplifier that attaches to any computer with a USB port. NuAmps can be used with Windows 98, 2000, XP and ME. NuAmps cannot be used with NT, because NT does not support the USB port. Beginning with SCAN 4.3, only Windows 2000 and XP are supported.

- Analog inputs: 40 unipolar (bipolar derivations can be computed) per headbox
- Sampling frequencies: 125, 250, 500, 1000 Hz per channel for all channels
- Sampling method: all channels sampled simultaneously
- A/D resolution: 22 bits
- Full scale input range: ± 130 mV
- Input impedance: Not less than 80 MOhm
- CMRR = 100 dB at 60 Hz
- Input noise: 1.0 μ V RMS (6 μ V peak-to-peak)
- Bandwidth: 3dB down from DC to 262 Hz, dependent upon sampling frequency selected
- Interface Universal Serial Bus (USB), full support for Plug-and-Play technology
- Supported electrodes (as long as they do not exceed the dynamic range of the amplifiers, which is approximately ± 100 mVs):
 - Gold, Ag/AgCl, Carbon electrodes with Touch Proof (DIN 42-802) style connectors
 - QuikCap Ag/AgCl electrodes with Plastic DSUB37F
 - Tin electrodes (within acceptable saturation range)
 - Sintered silver electrodes
- Digital (TTL) inputs/outputs; 14 inputs and 2 outputs
- Quality control of electrode application

- Measurement of contact impedance (at frequency 30 Hz) in impedance mode
- Constant monitoring of connection during recording
- Isolation: Optical Signal Isolation
- Display: 16-letter LCD with background light, displaying amplifier status or electrodes with impedance greater than specified cutoff.
- Power supply and energy consumption: From USB (5V), in active mode current ≤ 500 mA, in standby mode current ≤ 20 mA
- Electric safety level: According to EN60601-1(type BF), IEC601-1
- Multiple NuAmps can be cascaded to acquire data from up to 157 channels

Cascaded NuAmps Systems

SCAN4.3.1 software supports synchronous recording from several NuAmps headboxes, including acquisition and impedance measurement. Such systems, combining two, three, or four individual amplifiers, are referred to as “Cascaded NuAmps Systems”. As with the single NuAmps system, the cascaded systems are designed for acquisition of EEG and long/middle-latency EPs in the research environment. Because of the additional hardware that is required, there are some important considerations with cascaded systems. These will be detailed throughout this manual. It is important that you note these differences.

Installing NuAmps

Installation of NuAmps software and hardware is very easy. The basic steps are:

- 1) Unpack the system components,
- 2) Connect the isolation transformer (if supplied),
- 3) Install the SCAN software,
- 4) Install the NuAmps hardware,
- 5) Configure SCAN for NuAmps, and

- 6) Perform a quick system test.


NOTE: The power supplies MUST be connected to the isolation transformer, if one was sent to you, and that must be plugged into a grounded wall outlet (3-prong). The ground should be a genuine earth ground, and it should be inspected by an electrical engineer.

1. Unpacking the NuAmps

The NuAmps has been shipped in a container designed to reduce damage due to shipping. Please retain this box and its contents in case you need to return the system for any reason.

Open the box and check for the following contents:

1. NuAmps amplifiers
2. 3' or 15' USB cable to connect NuAmps to the computer. (The A-type connector goes to the computer; the B-type of connector goes to NuAmps).
3. SCAN 4.2 (or newer) CD, if purchased (CD and software key may be shipped separately)

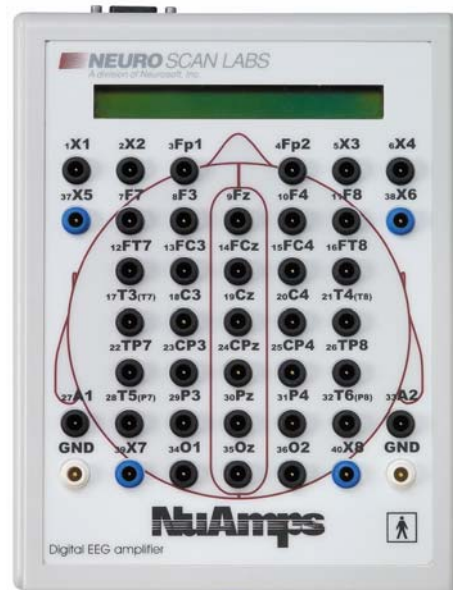

4. Software key, if the SCAN software was purchased (yours may be smaller than the one shown)



5. PocketTrace-2 unit and adapter



6. Adapter cable for the STIM-to-SCAN trigger cable. The 9-pin connector (DSUB9M) connects to NuAmps; the 25-pin connector (DSUB25M) connects to the STIM-to-SCAN cable.



Additional options include:

7. Isolation transformer, optional (separate box). If you purchased the laptop computer from Neuroscan, or if you run on batteries, no isolation transformer is needed. Some brands of laptops, especially those provided by Neuroscan, are powered by the portable power package. These computers can be purchased directly from other vendors, if desired (contact Neuroscan Technical Support for a list of these makers and models). All other configurations require the isolation transformer.

Cascaded NuAmps Systems

The following additional system components are included:

1. *USB-hub*. Several NuAmps amplifiers (up to 4) can be connected by a USB-hub. The physiological signals are synchronously registered from all amplifiers. Each amplifier comes with a USB-cable, PocketTrace, and a NuAmps-to-STIM cable. Individual USB-cables are used to connect the amplifiers to the USB-hub. (The hubs vary from vendor to vendor, and yours may not look exactly like the one shown below).



2. *Power Supply.* A Medical grade power supply with 5VDC, 5A output is provided. It powers the USB-hub.

2. *Connecting the Isolation Transformer*

Connect the computer power supply/transformer unit to the back of the computer. Connect the other end to the isolation transformer. If you purchased the laptop computer from Neuroscan, or if you run on batteries, no isolation transformer is needed. Some additional brands of laptops may not require the transformer (contact Neuroscan for details). All other configurations require the isolation transformer. (See the information below for Cascaded NuAmps Systems).



Warning

ALL components that connect to the computer should be powered through the isolation transformer.

3. *Install the SCAN Software*


SCAN software is required with NuAmps; NuAmps cannot be used with other software. All files needed to run NuAmps are included on the SCAN CD, and are installed automatically. The details for installing the SCAN software are presented in the Installation manual. Briefly, the steps for NuAmps are:

- 1) Connect the hardware lock, or dongle, into the parallel port on the back of the computer. This is the security lock - the software will not run if this is not in place. **The dongle is essentially your license, and it is the most valuable part of the system. Guard it carefully. The replacement cost is the cost of the complete system.**
- 2) Start you computer. If you have Windows 2000, *you must install and run the software when logged on as the administrator.*

- 3) Install the SCAN software from the CD. The installation process should begin automatically after you place the CD in the CD drive. If it does not start automatically, run "d:\setup.exe" (or whatever your CD drive is).
- 4) Reboot the computer at the end of the software installation, when prompted.

4. Install the NuAmps Hardware

Locate the NuAmps unit, the adaptor for the STIM-to-SCAN cable (if you have STIM), the USB cable, PocketTrace and its adapter, and the amplifier stand.

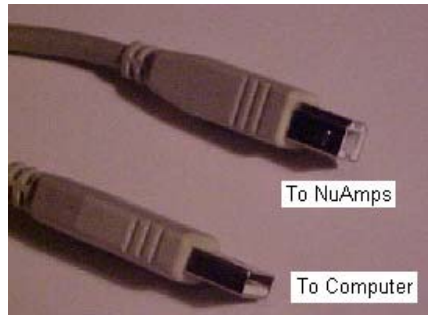

1. Connect one end (the B-type connector) of the USB cable to the NuAmps unit, and the other end (the A-type connector) to the USB port on the host computer.

2. Connect the PocketTrace signal generator to the Cap Connector on the amplifiers, if desired. The large gray connector plugs into the cap connector, and the longer black cable plugs into the jack on the side of the PocketTrace unit (use the mini-phone jack connector that fits the cable plug). Disregard the loose touch-proof cable connectors for now.

PocketTrace may be used in a variety of ways, as explained

more fully in Appendix C at the end of this manual.

3. If you have a STIM system, you will have a STIM-to-SCAN cable. This sends TTL trigger pulses from STIM to SCAN. With NuAmps, there is a short adapter cable that connects the SCAN end of the STIM-to-SCAN cable to the serial connector (DSUB9) on the top of NuAmps.



Warning


Be sure to connect this cable to the trigger port on the NuAmps unit; do NOT connect the STIM-to-SCAN cable directly to the parallel port on the SCAN PC. Secure the DSUB9 connector with screws to avoid a loose connection. Do not connect any other device to the serial port except for the STIM-to-SCAN cable, as other devices could damage the port.

Cascaded NuAmps Systems

You must install the NuAmps amplifiers one at a time, specifying the driver each time. Do NOT connect all the amplifiers to the hub, then connect the hub to the PC and turn it on. If you do, only the first NuAmps will be recognized. Please perform the following steps to set up a cascaded NuAmps system:

1. Connect the USB-hub using the 15” cable to PC. Then, connect power to the USB-hub with the provided medical grade 5VDC, 5A power supply (from AULT). Windows should recognize the hub and install the drivers automatically. Some hubs come with a CD that includes the drivers. Use that CD, if your computer requests you to do so. (Windows XP usually does not require that; Windows 2000 may make such a request).

2. Connect the first NuAmps to the USB-hub. Most USB-hubs come with 4 ports that are marked by

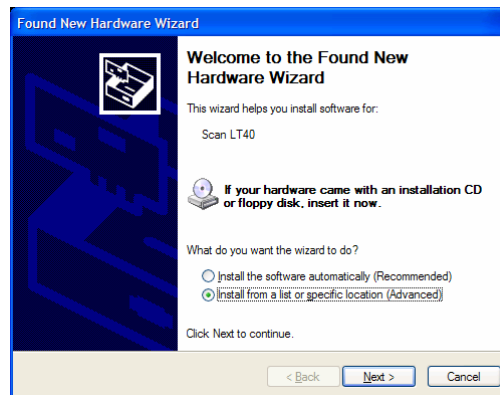

numbers from 1 to 4. Use the first port.

3. Windows will acknowledge that new hardware “SCAN LT40” is found, and try to install the driver automatically for a few seconds.

4. Because you are connecting the amplifier through the USB-hub, the automatic driver installation is not possible. As a result, the “Found New Hardware Wizard” dialog will appear.

5. Select “Install from a list or specific location (Advanced)”, and click the “Next” button. On the next dialog screen, select “Don’t search. I will choose the driver to install” radio-button, and click the “Next” button again. The OS shows the list of compatible hardware. Select “Scan LT40” and click the “Next” button. If the installation requests the “MCSUSB06.SYS” file, browse to your “Windows\System32\Drivers” directory (where the file was placed during installation), and click the OK button to proceed. The driver will be installed.

6. Connect the second amplifier to second USB port. The PC will again ask for the device driver. Repeat the sequence of operations described in steps 4 & 5.

7. Install the other amplifiers in a similar fashion. This may seem tedious, but once you go through the installation, Windows will recognize and load the drivers automatically thereafter.

NOTES

1. Each USB-connector on your PC has a unique number that is programmatically added to a newly attached device. Therefore, if you connect the USB-hub to another USB-port on your computer, the internally generated IDs for that amplifier will be different and you may need to repeat steps 4-7.


signals. Basically, the number of channels is a multiple of the number of headboxes and 40, minus the number of links between reference pins, which in turn is equal to the number of amplifiers minus 1. As a result, we get the following formula:

$$(N \text{ of channels}) = (N \text{ of amplifiers}) \times 40 - (N \text{ of amplifiers} - 1)$$

and table:

Number of cascaded units	Number of independent channels
2	79
3	118
4	157



5. Configure SCAN for NuAmps

When the computer has booted into Windows, the next step is to install the amplifiers using the AmpInst.exe program. The program is executed by any of the conventional methods in Windows. For example, go into the Windows Explorer, then to the Scan 4.3 folder, and double-click the  AmpInst.exe line. For future



convenience, you may want to create a shortcut to the program on your desktop: Amplifier Install . The Scan 4.3 AmpInst display will appear.

Select NuAmps and enable any of the Video Camera option if you have one installed and plan to use it with the SCAN software. Click OK. You then see a confirmation screen showing your selection. Go Back if need be, or click Finished.

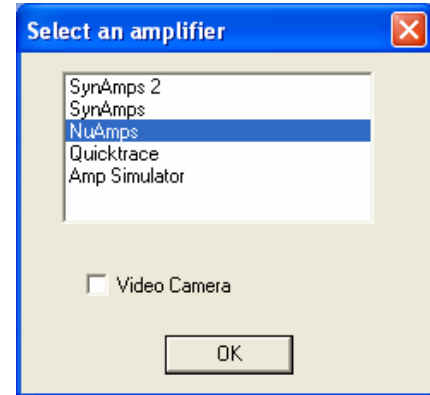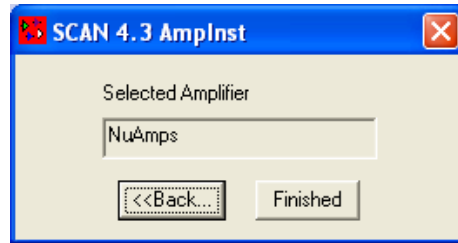

This step writes configuration information to an initialization file (you will not see anything happen).

Cascaded NuAmps Systems

The only difference in configuring SCAN for a cascaded NuAmps system is the setup file you select. In the ACQUIRE program, click “File->Load Setup...”, and select an appropriate setup from “Setup Files\NuAmps” subdirectory. There are setup files for 80 and 120-channel recordings, which correspond to acquisition from the specially designed QuikCaps. The recording will be performed from 2 or 3 amplifiers, respectively. SCAN 4.3 contains 2D MAP files for cascaded systems.

If you need to modify the setup file or create your own, then follow general instruction of setup file creation and modification. Remember about the number of independent channels, and to create the appropriate 2D map file in which the physical channel numbers and output labels are matched (as described in the ACQUIRE and MAPGEN manuals).

Differences Between Single and Cascaded Systems


Please be aware of the following differences between single and multiple NuAmps systems:

1. Channel usage and numeration are consecutive. Channels from the first amplifier are used first, then from second amplifier, third amplifier and, finally, the fourth unit. For example, consider the case of a 50-channel acquisition. SCAN will acquire data from 40 channels from the first amplifier, numbered from 1 to 40, and then acquire 10 extra channels from second amplifier, numbered from 1 to 10. There is no way that you can record channels in any other order (all of the channels in one NuAmps MUST be used before adding channels from the next one).
2. If you activate the “NuAmps Setup” dialog screen, the “Amplifiers” page will show you the sequence of amplifiers.

This sequence coincides with the port numbers from the USB-hub to which the amplifiers are attached. This sequence specifies which amplifier will be used first, which is second, etc.


The following information is displayed for every amplifier:

- Model number with firmware version (latter should be at least 4.4);
- Number of channels (always 40 for NuAmps);
- Serial number (as indicated on the back of the amplifier);
- Hardware version (currently 3 for the NuAmps DC amplifier).

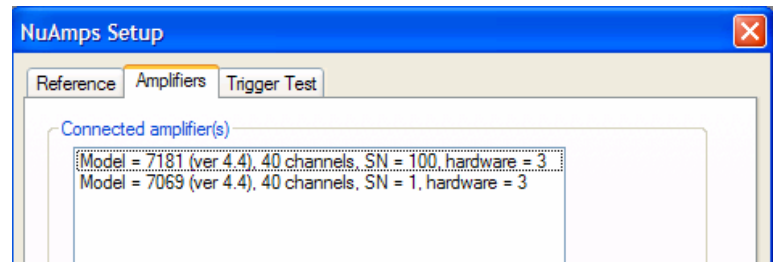


The serial number in this list can help you identify which amplifier is first in chain. Remember that you **MUST** connect the NuAmps-to-STIM cable to the first amplifier to be able to record TTL trigger impulses. If you connect the trigger cable to any other amplifier, then no trigger information will be available.

3. As you start recording, the NuAmps LCD display also shows the ordinal number of the amplifier. It is an additional reminder for you to identify amplifiers in order, and to attach the NuAmps-to-STIM cable to the proper amplifier.
4. Select a setup file with more than 40 channels, then activate “NuAmps setup” and go to the “Reference” tab. You will see the following warning at the bottom:

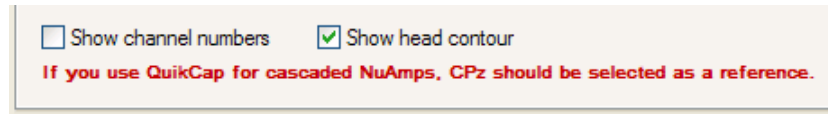

The CPz electrode is wired in the Quik-Cap for use as the Reference electrode across headboxes. The GND electrode is also "hard-wired" in the cap to go to all headboxes. This ensures that potentials recorded by different amplifiers are in the same “coordinate” space. Therefore, different analysis, like PCA, ICA, or source localization, can be applied validly to the recordings.

5. If you are using individual leads, you can use any Reference, but remember to link together the corresponding Reference pins on all amplifiers. Do the same for the GND pins.

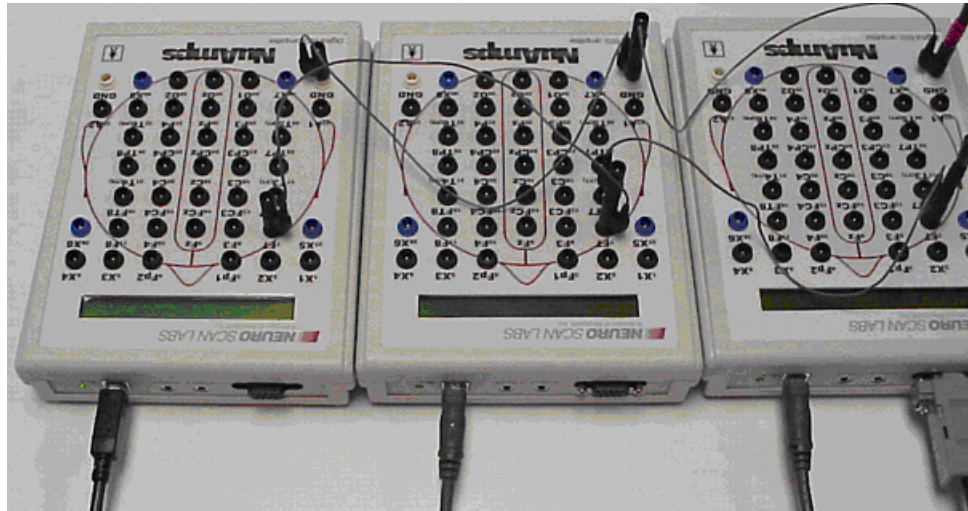

6. In the figure above, the F7 electrode was selected as the Reference. You can see that all F7-pins are chained, as are the GND-pins (because there is pair of GND pins per amplifier, you can use either of them - they are linked physically inside the headbox). The amplifier at the right side is the first amplifier in the chain, and therefore the NuAmps-to-STIM cable is connected to it.


Acquisition from Cascaded NuAmps

Acquisition and impedance measurement for cascaded NuAmps are quite similar to acquisition from a single unit system, with several important differences:

1. Before acquisition begins, the ACQUIRE programs the amplifiers to work synchronously, and that usually requires several seconds. During that operation, ACQUIRE shows the “Programming amplifiers” message.
2. If you select a setup file where the number of recorded channels exceeds the number of physical channels that are present, the program warns:

To resolve the problem you need either to decrease number of channels or attach an additional amplifier.

3. Supplied setup (AST) files are configured in such way that you can use select different display pages to see signals from each amplifiers separately. This is especially useful when recording from three or four amplifiers.

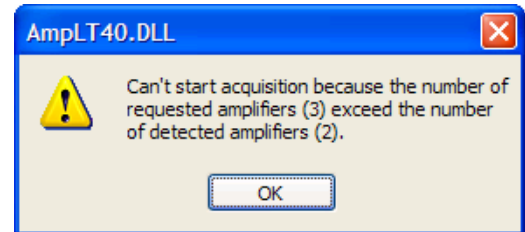

Tuning Performance of Windows XP Systems and Cascaded NuAmps

Recording from several NuAmps units increases the PC load. If you are using a slow computers under Windows XP, especially with the initial release of XP, then recording from cascaded units may overload your system. Overloaded systems become less responsive to your actions, have slower responses and may even crash. This is especially likely if you are using three or four amplifiers on a slower computer (<1.6 GHz), with the original release of XP installed.


To check for that situation, please do the following:

1. Right-click on the taskbar and select “Task Manager” from the popup-menu.
2. Switch to the “Performance” page and verify that “Show Kernel Time” is checked under the “View” menu.
3. Now start acquisition from the cascaded system and see if the green line, showing the overall system load, reaches 100% and stays there for several seconds, minutes or until you stop acquisition. Or, your system may show a huge load, always at level 90% as on the picture below:

To tune your computer, you need to install updates from Microsoft; especially the ones called “Service Packs”, or those related to the USB. Windows XP made this updating easy. Click on the “Start” button on your taskbar, select “Help and Support”, and then click “Windows Update”. Your computer will connect to Microsoft web-server to get the latest updates. Follow the instructions.

After updating, test the performance on your system again. If you still get an overload, you might need to upgrade your hardware. Try increasing the RAM (up to, for example, 256MBs), and obtain a faster graphics card. For graphics cards, try enabling full hardware acceleration to improve the overall performance. If your system still is overloaded while acquiring from the cascaded system, you should get a faster computer, especially one with a CPU faster than 1.6 GHz.

Interestingly, Windows 2000 systems are less demanding, and computers with CPU speed such as 900 MHz show the load level at 30% for 4 cascaded NuAmps!

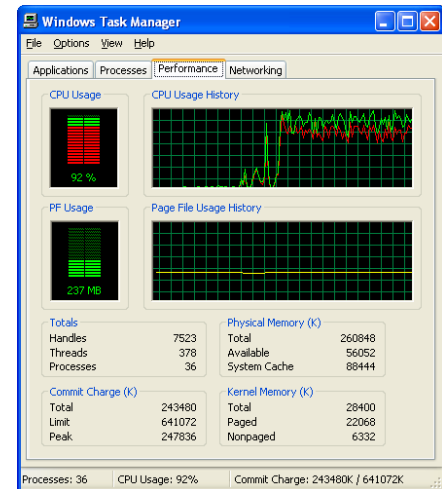

Safety and Quality Warnings


1. A Medical grade power supply **MUST** be used to power the USB-hub to ensure:
 - Patient safety.
 - Sufficient current for 4 amplifiers. (Because each NuAmps requires 0.5A current, then 4 units require 2 A. Some standard USB-hub power supplies cannot provide 2A current, even though they may claim to do so).
 - Reduced instrumental noise.
2. Isolation transformer **MUST** be used for desktop systems.
3. Portable power package **MUST** be used for portable systems.
4. Plug your STIM (or compatible) system to the same wall power outlet.
5. The same Ground and Reference electrodes should be used across amplifiers. Quik-Caps are hard-wired in that way, so this is accomplished for you. If you are using individual electrodes, then you must link the Ground pins across all amplifiers and the Reference pins across all amplifiers.
6. Trigger information is recorded from the first amplifier in cascaded system; the NuAmps-to-STIM trigger must be connected to the first NuAmps in the chain.




Failure to do so may affect patient safety (per warnings 1-3), reduce the quality of recorded signals (per warnings 4-5), or block some operation (per warning 6).

6. *Quick System Test*

Click the Scan 4.3 icon on the desktop , and see the Program Launcher.

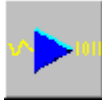

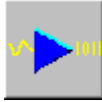

Click the ACQUIRE icon , and the main screen in ACQUIRE will appear. Click the Open File icon on the Toolbar , and select the NuAmps40.ast setup file (from the \Scan4.3\Setup Files folder). Then click on the green arrow on the Toolbar  to begin acquisition. The signals from NuAmps will appear in ACQUIRE (flat lines are OK), and you will see "NuAmps start" on the LCD display. If the LCD does not light up, or if you get any error messages, please contact Technical Support.


The Amplifiers

Each NuAmps amplifier provides inputs for up to 40 monopolar channels, per headbox. All channels are always recorded as true monopolar channels. The active leads are referenced to the ground lead, including the reference lead(s). The reference(s) you select (from the Reference tab on the NuAmps Setup dialog box) is applied in the software. Bipolar channels are derived mathematically from the monopolar channels (described below).

You will see in the face of NuAmps the standard electrode placements, with some auxiliary ones, totaling 40 inputs.

The connector for the electrode cap is on the lower side. There are several connectors on the upper side of the unit. In addition to the USB connector, mentioned above, there are the SYNC1 and SYNC2 jacks, and a 9-pin DSUB9F connector (for use with the STIM-to-SCAN cable).

SYNC1 and SYNC2 are alternate trigger inputs for single bit TTL pulses. These are useful in configurations that include simple stimulus presentation systems other than STIM, that have one or two bit outputs. A TTL pulse in SYNC1 will result in a trigger type code of 1 in the CNT file; a pulse to SYNC2 results in a 2 in the CNT file. The TTL pulses should be 5V, positive logic, with a duration of approximately 5-10ms. The jacks are 3.5mm stereo connectors with 3 contacts: 1) TTL Out, 2) TTL In, and 3) USB Ground. The TTL inputs/outputs have galvanic isolation from the patient.

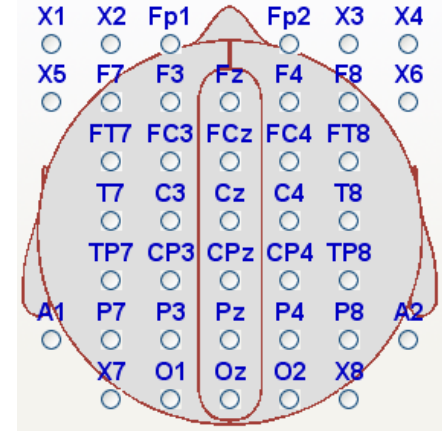

Software Control of NuAmps

Data acquisition on the NuAmps is controlled by the host computer executing the ACQUIRE program. Acquisition, control, display and storage of data are all performed by this program. For a detailed explanation of this program, please

see the ACQUIRE section of the SCAN manual. The sections below describes those aspects of ACQUIRE that are specific to NuAmps.

NuAmps Setup Files

To record signals from NuAmps, you will need to use one of the supplied setup files, modify one of the supplied ones to fit your needs, or create a new setup file from scratch. Setup files contain all the system configuration information, and are stored as separate files (with an AST extension). The following NuAmps setup files are installed in the \Scan4.3\Setup Files folder:

NuAmps 19.AST - 19 standard electrodes from 10-20 system:
(Fp1,Fp2,F3,F4,C3,C4,P3,P4,O1,O2,F7,F8,T3,T4,T5,T6,Fz,Cz,Pz).

NuAmps 30.AST - 19 standard electrodes as in NuAmps19.AST plus 11 additional electrodes: FT7, FC3, FC4, FC4, FT8, TP7, CP3,CPz,CP4,TP8, Oz. Notice that some electrodes from the 10-20 system are labeled according to the extended 10-20 system, such as, T3->T7, T4->T8, T5->P7, T6->P8.

NuAmps4 0.AST - includes all 40 electrodes found on the amplifier with labels that correspond to the amplifier overlay.

NuAmps 40 QuikCap40.AST - setup file to be used with NuAmps QuikCap. Basically this is the same as NuAmps40.AST, except 4 electrodes are renamed according to the cap: X5->FT9, X6->FT10, X7->PO1, X8->PO2.

NuAmps 80.ast - Setup file to be used with cascaded NuAmps system.

Most of the setup files have the following parameters:

Sampling rate - 500 Hz

Filters: HPF - 0.5 Hz, LPF - 70 Hz, Notch - off

Reference - A2.

Much of the information in the ACQUIRE manual is relevant for creating setup files with NuAmps, and most of it will not be repeated here. The unique displays seen only for NuAmps are described below.

Whenever possible, you should use one of the supplied setup files; however, in some circumstances you may need to create your own. To demonstrate the process, we will create a 36 channel setup file from scratch. In ACQUIRE, select Edit (on the Main Menu bar), then select Make Default Setup, and click OK to the warning. This enters the default settings throughout the setup file.

Enter Number of Channels. Click Edit, then Overall Parameters, then the Amplifiers tab. (You may see two messages regarding high and low pass filter adjustments). Begin by entering the number of

channels you wish to record , and then click the Reset Positions tab.

This positions the electrodes in a 6 x 6 matrix seen on the right side of the display.

Note that there is only one gain setting for NuAmps (19) which corresponds to an input range of ± 131.5 mV (and 22-bits). The resolution is $131.5/(2^{21})=0.063\mu\text{V}$.

Gain	19	▼
Range	263 mV	
Accuracy	0.063 $\mu\text{V}/\text{LSB}$	


1	2	3	4	5	6
7	8	9	10	11	12
13	14	15	16	17	18
19	20	21	22	23	24
25	26	27	28	29	30
31	32	33	34	35	36

Range and *Accuracy* information are related to the Gain setting (there is only the one option with NuAmps).

The Range values indicate the voltage limit ($\pm 131.5\text{mVs}$, or 263mVs , peak-to-peak) for each channel - if the incoming voltage exceeds these limits, the channel will saturate. Accuracy refers to the precision of voltage measurement along the y-axis. The value displayed indicates in microvolts the resolution ability. If the accuracy is, for example, $.063\text{ uV}$, then the voltage resolution will be in $.063\text{ uV}$ steps. This is similar to the dwell time on the x-axis, or, in other words, the *time* difference between adjacent data points. Accuracy is the least measurable *voltage* difference between points.


All of the other settings on the Amplifiers screen (and for the Overall Parameters settings) are as described in the ACQUIRE manual, with these exceptions.

AD Rate. There are 4 AD rates available for NuAmps: 125, 250, 500 and 1000Hz

A/D Rate . The 125Hz rate is not available with Cascaded NuAmps Systems.

AC/DC. The only option for NuAmps is DC.

Filters. The Low Pass filters do not extend as high as for SynAmps (due to the maximum AD rate of 1000). The available filter options are dependent on the sampling rate, and cannot exceed the Nyquist frequency (half of the sampling rate). As a result, if you decrease the sampling rate, the program may change the low pass filter automatically to the default 70 Hz value. Make sure this still fits your needs. High pass values depend on the sampling rate.

32-Bit Acquisition. The  32-bit acquisition option lets you acquire data in either 32-bit or the original 16-bit modes. If you are recording in DC, you must use 32-bit acquisition because of the wider dynamic range it allows, and to record the slowest frequency activity accurately. In general, 32 bit acquisition is preferred in all cases. The one consideration is that it will essentially double the file size as compared to the 16-bit mode.

Triggers tab. With NuAmps, the only available triggering option on the Triggers tab is External (no Internal or Voltage triggering).

Enter Channel Labels. When you create a setup file from scratch (such as, after selecting the Make Default Setup option under Edit), you will see numbers instead of the common channel labels. To change these, click Edit, then the Channel Assignment Table (shown in part).

Assignments				Physical Channels			
	#	Phys	Label				
+	1	1	1	1	2	3	4
+	2	2	2	9	10	11	12
+	3	3	3				
+	4	4	4				

Here we need to decide the order of the channels as they will appear on the display. Let's say we want the first channel to be FP1. Note on the face of the NuAmps amplifiers that FP1 uses the #3 input, or third physical amplifier channel. Therefore, enter 3 in the Phys. Chan. field, and FP1 for the label.

Continue in this way to enter all of the electrodes and their respective physical channel numbers, referring to the amplifier face for the numbers. Use the electrode labels shown on the face of the amplifiers, insofar as that is convenient. You can use different labels; however, there are other sections in the program that read the labels for automatic positioning (and other operations, including mapping and source localization in SOURCE). If you use nonstandard labels, these will not be recognized.

	#	Phys	Label
+	1	3	FP1

The final Channel Assignment Table appears as follows (shown in part).

In this example, we are recording A1 and A2 as regular channels, and X1 and X2 are VEOG+ and VEOG-, respectively, and X3 and X4 are HEOGL and HEOGR, respectively (not shown above). When you are through, click OK (and save the setup file, if desired).


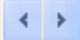
Channel Layout. Under Edit, select the Channel Layout option (the complete description of the Channel Layout dialog box is found in the ACQUIRE manual). Channel Layout is used to position and size the electrode displays. In this example, the positions are in a 6 x 6 matrix. You can reposition the displays manually if you want, but it is easier to do it automatically (see the ACQUIRE manual for complete details).

Assignments			
	#	Phys	Label
+	2	4	FP2
+	3	7	F7
+	4	11	F8
+	5	8	F3
+	6	10	F4
+	7	12	FT7
+	8	16	FT8
+	9	13	FC3
+	10	15	FC4


For now, retrieve the setup file, and go to Channel Layout. Click the **10-20 Positions** button and then the **Match Labels** button. This recognizes the standard 10-20 system labels, and positions the electrodes accordingly. If you used different labels, will need to position those channels manually. Then click OK. The new positions will be transferred to the channel displays.

For further adjustments, click the **Adjust Positions** button, and see the following screen.

The three groups of adjustments are for Spacing, Size, and Translation. **Spacing** expands or compresses the grouping of the displays either


vertically , or horizontally . Similarly, the **Size** buttons increase or decrease the size of the electrode displays either vertically or horizontally. The **Translation** buttons will shift the entire grouping of the displays vertically or horizontally.

When you have the displays positioned and sized as desired, click OK, and then OK again. (It is a good idea to re-Save the setup file to save any changes).

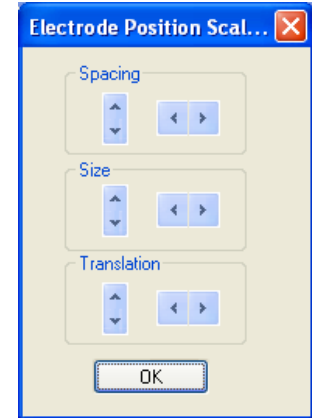

NuAmps Setup Screen. Click Edit (on the Main Menu bar), and select the "NuAmps Setup..." option. You will see the following display, with three tabs.

Reference. The Reference display lets you select what site(s) to use for the reference electrode. The NuAmps amplifiers are different from most other amplifiers. The main point to understand is that the reference is computed in the software. That means there are few limits on what can be used for the reference. You can select a single electrode, used linked or single ears, and so forth. Let's say you want a linked ears reference. That is the first option (A1+A2)/2. Click that button, and you will see the current selection written at the top of the display.

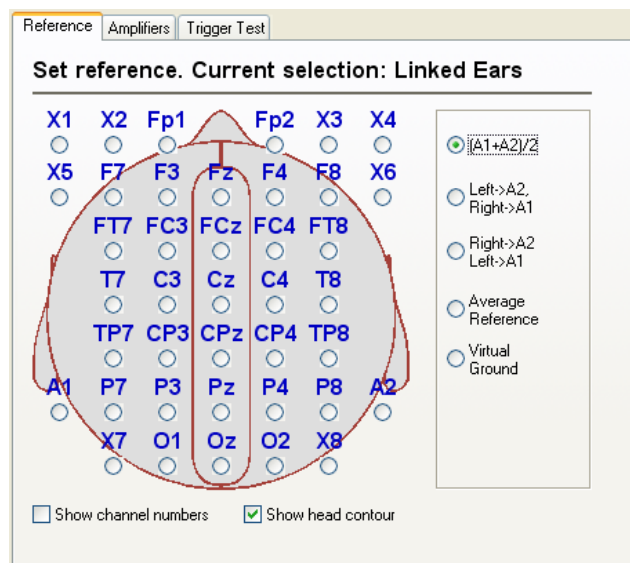

Set reference. Current selection: Linked Ears

As mentioned, you have a great deal of flexibility in the selection of the site(s) for use as the reference. To select a single channel, such as Cz, just click the Cz field. You will see the choice indicated in the top line of the display.

Set reference. Current selection: Cz

The complete list of reference options are:

(A1+A2)/2: This is the linked ears option.

Left=>A2, Right=>A1: This option uses the contralateral ear/mastoid as the reference.

Right=>A2, Left=>A1: This option uses the ipsilateral ear/mastoid as the reference.


Average Reference: In this case, all specified electrodes are recorded and the averaged potential is calculated and subtracted from original signals. Only actually placed electrodes are used. For example, if you specify 10 electrodes to be recorded, but only place 8 of them, then only this group of 8 electrodes will be used.

Virtual Ground: This option should be used only for demonstration purposes while recording signals from the PocketTrace signal generator. In that case, the program displays pure potentials under active electrodes relative to ground. Such traces, however, may have 50/60 Hz noise, DC drift, etc. The Virtual Ground can be used only for high amplitude signals such as ECG. For low amplitude signals, like EEG, differential montages are used, and the reference potentials (for example, the earlobe) are subtracted from the active potentials.

Note that you have the options to display the head contour and the channel numbers at the bottom of the display.

Amplifiers. Click this tab to query the NuAmps headboxes you have attached.

If you are recording in the 16-bit acquisition mode in an fMRI environment, it is recommended that you enable the "Use moving average" option, and then select a window from 1-10 seconds. The option is not needed with 32-bit acquisition.

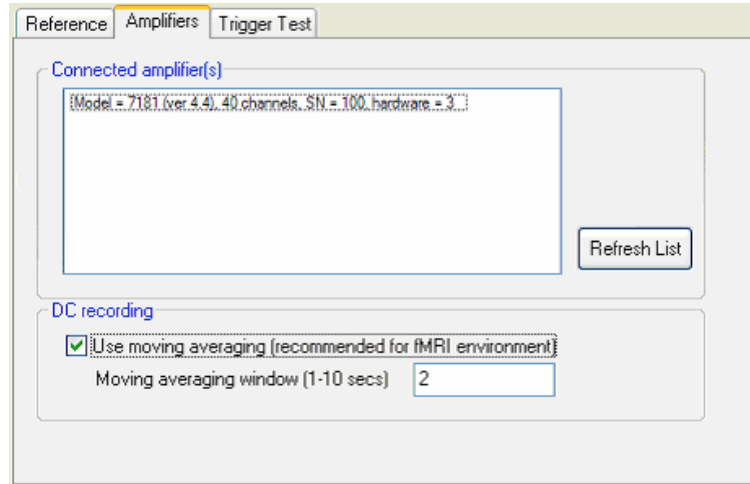

Trigger Test. The Trigger Test option is used for diagnostic purposes if you are having triggering problems. It reads the status of the serial port on the NuAmps trigger input. Its operation is described in Appendix A at the end of this manual.

This completes the basic steps for creating a setup file for NuAmps (save the setup file). You might also want to create bipolar and other re-montaging files. This process is described next. If you purchased a photic stimulator, you will also need to set up the recording sequence. This is described in the Photic Stimulator manual.


	Decimal	Binary
Stimulus	0	00000000
Last Stim	0	00000000
Response	0	00000000
Counter	0	
Hold Value	0	

Test


Creating a Bipolar Montage. NuAmps records all channels as monopolar channels, referenced to ground. The Reference you select is applied in the software. Bipolar montages are computed from the referenced EEG signals by simple subtraction. The Montage Editor is used to create an LDR file (linear derivation file) that performs the subtraction. The complete operation of the Montage Editor is explained in an appendix to the EDIT manual. We will demonstrate here how to create a simple LDR file that can be used to compute selected bipolar channels. (You can create, modify, or apply the LDR file online in ACQUIRE or offline in EDIT. We will demonstrate the online option here).


In this example, we have a VEOG+ electrode (above the eye) and a VEOG- electrode (below the eye), and we want to create a bipolar channel: VEOG+ minus VEOG-.

1. In ACQUIRE, retrieve the setup file you will be using.

2. Click the Montage Editor (ME) icon from the Toolbar .

3. In the Montage Editor, click File, then Create New Montage. Use the same number of channels, but make sure you click the Create Identity Matrix field. This creates an output LDR file that is identical to the original one. The matrix will show a 1 in each cell as you look down the diagonal from upper left to lower right cells. Click OK.

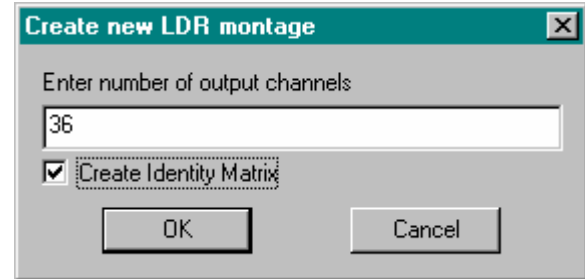

4. Unclick the Show Head Contour icon . This lets you see more of the matrix.

5. Right click on the VEOG+ label in the column on the left, and select Rename channel.. Enter BipVEOG (for bipolar VEOG). *Do NOT use spaces in the new label names.*



6. Right click on the VEOG- label, select Delete channel, and delete the VEOG- channel (ignore the warning message). This is not mandatory. We are setting up the LDR file to show the single bipolar artifact channel. We could also keep the existing VEOG channels and create an additional bipolar one.

7. Now, look on the row corresponding to the BipVEOG channel we are creating. Look at the intersection of the BipVEOG row, and the VEOG+ column. There should be a 1 in that cell already (created when we selected

the identity matrix option). Now find the intersecting cell for the BipVEOG row and the VEOG- column. *Add a -1 in the cell in the VEOG- column.*

This says to take each original data point from the VEOG+ column and multiply it by 1 (no change), and to take each data point in the VEOG- column and multiply it by -1. The new channel, BipVEOG, will be the linear sum of the data points for all channels in the row. Each new data point will be VEOG+ minus VEOG-, or the derived bipolar channel. (The other channels are all zeroes, and do not contribute).

(0:0)	F7	FP1	VEOG+	VEOG-
F3	0	0	0	0
FP2	0	0	0	0
F7	1	0	0	0
FP1	0	1	0	0
BipVEOG	0	0	1	-1


8. Click File, then Save Montage File. Enter a file name, and *be sure the "Save as type:" field is set for Linear Derivation Files.*

File name:	BipVEOG	Save
Save as type:	Linear Derivation Files (*.ldr)	Cancel

The LDR file that is created can be applied anytime you use the same setup file, or with any data file recorded with that setup file. Once the LDR is created and saved, you do not have to create it again. Close the Montage Editor.

9. In ACQUIRE, when you are acquiring data (Single Window Display only), and wish to see the derived bipolar channel, all you need to do is right click in the data display and select Add Derivation.

Select the LDR file you just created, and you will see that data displayed with the new channel.

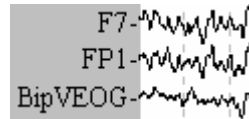

Click Remove Derivation (right click menu) to restore the data to its original form. In the offline EDIT program, you can apply the LDR file in the same way, or you can apply it with the Linear Derivation transform (to save the new file).

DC Electrode Considerations

The DC capability of the NuAmps amplifier presents special requirements that most electrophysiologists more familiar with AC coupled systems may not have encountered. With DC amplifiers, battery potentials generated from electrodes are not dissipated by the decoupling capacitors found in AC only systems. NuAmps cannot correct offset voltages greater than plus or minus 130mVs. Electrode combinations generating battery potentials greater than 130mV will saturate the amplifier. Saturation will become evident when the amplifier displays a flat line with no apparent activity and a saturation value greater than 100% is displayed. Under normal operating conditions you should not experience saturation values greater than about 30%. **If your electrodes show a consistent DC shift that climbs quickly to saturation, then your electrodes are not suitable for DC recordings.** There are several rules to follow when selecting electrodes to avoid these problems:

- 1. Never use different combinations of paste/gel and or electrodes.** For example, a common error is to use one type of electrode and paste in a monopolar derivation (i.e., electrode cap and gel) and a

different electrode and or paste on the reference lead. By using different metals and electrolytes a battery potential has been created and the electrodes will drift. If you have the same electrode materials and gel and you are still experiencing saturation, then double check your electrodes. They may not be made out of identical material. Measure the DC voltages across these leads and you will possibly find a large offset potential.

- 2. Keep the interface between the skin and electrode consistent.** For example, with an electrode cap the interface with the skin is a nonconductive rubber and the conductive gel. The gel makes contact with the electrode. If you place a reference electrode directly on the subject's skin for a reference, you have created a different interface. You now have a gel plus metal to skin interface. The best way to avoid this problem is to obtain a separate but identical electrode to employ as a reference. Another technique is to separate the metal from the skin with an adhesive electrode collar. *Note - drift problems with the reference electrodes in multichannel recordings are usually observed across all monopolar derived electrodes.*
- 3. Select metals that are known to produce the smallest battery potentials.** Sintered Ag/AgCl electrodes are probably the best and easiest to obtain. If you use these electrodes, remember to take care not to scratch the surface.
- 4. The DC level is relatively independent of the electrode impedance.** If you have set your electrode impedances to the standard 5kOhms or less and your electrode is still drifting, additional work on the impedance will usually have no effect (except on the subject). Since the primary source of battery potentials is the interface between the gel and electrode and not the skin, further work on the gel to skin interface will probably not help. This is a good time to replace the electrode or examine the metal to wire solder joint (another potentially large source of battery potentials).

5. Record in a comfortable and cool environment. Sweat potentials can be a major problem for DC recordings. They produce transient and unpredictable results.

On some AC coupled systems, high-pass values below .1 Hz are often available as an approximation to DC. If you are interested in slow potentials, there is no need for these values since the NuAmps is a DC coupled system. The advantages to recording slow potentials with a DC high-pass are two fold.

First, **DC amplifiers are less prone to impulse artifact.** The recovery time of an AC coupled system (actual or simulated as is the case on the NuAmps) can be described by the time constant of the system. The time constant is the time needed by an amplifier to decay to 37% of the peak response to a sudden impulse. To calculate the time constant of a simple filter use the following formula where T is the time constant and f is the cutoff frequency:

$$T = \frac{1}{2 \cdot \pi \cdot f}$$

This formula can be used to determine the effects of high-pass filter values you may have on your recordings. For example, if the amplifier encounters an impulse artifact with a 0.01Hz high-pass filter, it will require 16 seconds before the system returns to 37% of the original amplitude! If you are using such a filter setting so as to approximate DC, we strongly recommend that a DC value be employed. The DC setting is insensitive to impulse artifact. The NuAmps is not affected by even the most substantial of artifacts such as those generated by somatosensory and magnetic stimulators when recording in the DC mode.

AC coupled systems that employ long time constants such as .01Hz can be very unstable and are highly susceptible to sudden impulses (movement artifact) causing the subsequent smoothing of the impulse, consequently, these amplifiers saturate and need an external reset circuit to restore the system.

Second, a **DC recording can always be refiltered with different high-pass values.** Digital filtering can be applied to the data to examine the effects of different high-pass values. This is not true,

however, if the data have been acquired with an AC coupled system. You will have to accept whatever high-pass values were originally used to sample the data. In addition, digital filtering (in both EDIT and ACQUIRE) can offer sharper frequency cutoffs without the phase shift of an analog filter.

Notch Filters and Noise

Narrow band notch filters centered at 50 and 60 Hz are available to reduce main power frequency interference. Although these filters have a +/- 1Hz window, and affect a narrow range of the EEG spectrum, caution should be exercised in their use. Many evoked potentials have energy in the 50-60 Hz band and data may be significantly distorted. Notch filters should not be used routinely. Rather, attempts to reduce noise in the recording environment should be made first before adding a notch filter.

The option for the filter is found on the Amplifier screen under Edit, under Overall Parameters..., from the Acquire Main Menu bar.



60 Hz (50 Hz) Interference - Here are a few quick things to check if you have significant interference in your data:

Impedance of leads - Impedances of all leads should be below 5kOhms.

Impedance of ground lead - It is important that a solid ground lead be placed on the subject. The ground lead is used by the amplifiers to reject common-mode interference such as 50 and 60 Hz main frequencies. Double check this lead to be sure it also is under 5kOhms.



Warning

If all of the channels show high impedance values, this is an indication that the Ground or Reference electrodes have high impedances themselves.

Power cables - Make sure that there are no power or video cables near the electrode leads, the subject, or the amplifiers. Position the USB cable so that it does not cross the electrode leads. Although signals coming from the amplifiers are amplified, draping a power cable over the amplifiers or amplifiers cable will produce interference.

Power supply - The power supply for your laptop computer can introduce line noise. If you encounter excessive noise in the recordings, try running on battery alone. If that resolves the problem, it may be that the power supply is just noisy, or, there is not a true earth ground.

Monitors - Video monitors are now a common device in the neurophysiological laboratory. Unfortunately, they can also be a serious source of interference, especially the older ones. Most monitors radiate more from the sides than the front. If you suspect interference from your monitor, try moving the orientation of the screen. It should vary with orientation. Another clue that a monitor is emitting noise is a clear pulse with an interval corresponding to the refresh rate of the screen (i.e., 16.667 ms for a 60Hz refresh rate). It should be noted that the newer low radiation monitors have significantly reduced noise radiation.

Other noise sources - Watch out for anything in the surrounding area of your laboratory that can generate an electromagnetic field. Any device with a large electric motor (centrifuges, freezers, elevators, large fans) or transformers (X-ray machines) is a potential problem. Also,

check to make sure that you have a good and 'quiet' ground connection within your building. A high quality milligauss meter can be used to measure these fields.

Calibration


Calibration of NuAmps is performed only in the factory, and no amplifier adjustments are possible on site. Inter-channel variations of $\pm 1\%$ are expected, even with high precision elements.

It is possible to measure the variation across amplifiers using an external signal from, for example, PocketTrace. Set it to produce a sine wave at, say 10Hz (see the PocketTrace Appendix for details). PocketTrace will deliver the signal to all channels. Measure the peaks in ACQUIRE, or, save the file, perform an FFT, and compare the peaks at the frequency you selected.

Impedance

Electrode impedances can be checked at any time while the subject is connected to the amplifiers. The leads are measured against the reference; the GND lead is used for recording signals. Therefore, all measurements for these leads will be affected by impedance of the reference and GND. (See Appendix E for more details). Follow these steps to check electrode impedance:




Step 1 - Place leads on the subject with standard electrode application techniques. If you want to see the impedances as you place the electrodes, be sure to place the ground and reference ones first.

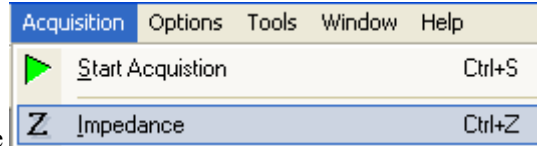

Step 2 - Click on the **Z** Impedance option under Acquisition from the Acquire Main Menu bar, or by clicking the Impedance icon **Z** from the Toolbar. The 'Impedance' display screen will appear.

The impedance of an electrode, or the opposition of AC current flow, is the result of the complex interaction between skin, electrolyte and electrode. Measured impedance will vary according to the transfer function of this junction.

Step 3 - Enter the threshold and maximum displayed impedance value. The default threshold value is 5kOhms, and the maximum value is set to 50 kOhms. To change these, click on the **Range...** button and enter the desired values. Then click OK.

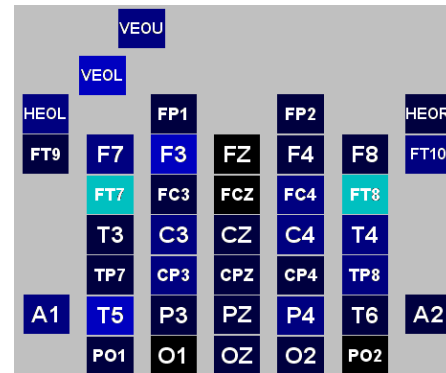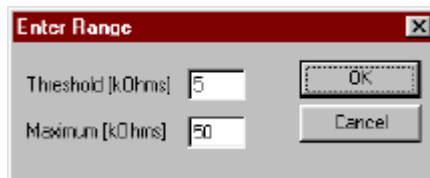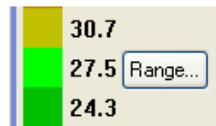

Step 4 - Impedance testing will begin when the Impedance screen is displayed. The impedance operation uses a 30Hz sine wave. The color-coded impedances will be displayed for each electrode according to the color of the display.

NOTE: The impedance of the Reference electrode cannot be measured by NuAmps, and it will always show pink in the impedance display, indicating maximum impedance.

Subject Ground

The subject ground should be placed as close as possible to the recording electrodes. Under most circumstances this location is somewhere on the head. The subject ground is used to reject common mode noise. Therefore, by placing the ground near the recording electrodes the common mode rejection ability of the amplifier is optimized. Avoid placing the ground at distal sites (arms and legs). An older technique often used to reduce stimulus artifact was to place the ground near the stimulating site. However, this procedure should not be used on NuAmps and will degrade system performance. Do NOT add a second subject ground. A ground loop could result, thus introducing 50 or 60Hz noise in the recordings.

Electrode Safety

The NuAmps amplifier inputs and attached electrodes are Type BF, which means in part that they are not connected to Earth Ground or Chassis Ground. Maintain this separation from Earth Ground by ensuring that the electrodes and any conductive parts of their connectors do not touch conductive parts, including the system enclosure or other grounded devices.

FAQs

The following are some of the common questions and answers regarding NuAmps.

1. Why I am getting requests to install the NuAmps driver for successive amplifiers?

The first time you install the multiple NuAmps, you must install them one at a time, and specify the driver for each when it is detected. See the cascaded installation directions above for details.

2. My computer reboots (or shows a blue screen) every time I start acquisition from the cascaded NuAmps system. What is wrong?

This is a symptom that your computer is being overloaded by the multiple amplifiers and cannot keep up with the demand. See the "Tuning Performance of Windows XP Systems" section above for suggestions. Rebooting can also happen if you do not have the power cable attached to the USB hub. In that case, the hub cannot pull enough current from the computer, and the computer reboots in self-defense.

3. My computer reboots if I disconnect NuAmps during acquisition?

USB devices in general, and the NuAmps in particular, should not be disconnected while in use.

4. Can I use a single PocketTrace to verify the operation of two or more amplifiers?

Yes. You can use PocketTrace individually for each amplifier, and verify the signals on the corresponding channels during acquisition. You can also jump the output from PocketTrace to additional headboxes. If you do that, you will need to connect the GND wire from PocketTrace to the second amplifier GND pin. Be sure to set the Reference to "Virtual Ground" whenever you use PocketTrace (and then set it back afterwards).

5. I have several USB-ports on my computer. Can I use them instead of the USB-hub and medical grade power supply?



This possibility is very attractive, especially for notebook systems with two NuAmps, as a way to reduce the number of components in the system. In some cases it is possible, depending solely on your computer. If you have only two cascaded NuAmps, try connecting both directly to the USB-ports of your computer. Then run ACQUIRE and start acquisition. If cascading is not possible, you will see the following warning:

If you did not see that message, and signals are scrolling across the screen, then you do not need the external USB-hub and its power supply.

To understand why it is possible in some cases and not others, you need to understand how PC vendors support multiple USB-ports on your computer. It is now very common to find two USB ports on laptops and four USB ports on desktop models. To provide communication through these ports, special USB chips are placed on the motherboard. If you open the “Device Manager”, and then expand “Universal Serial Bus controllers”, you might see something like the following:

In the above example, an IBM desktop system (with 4 USB ports) informs that it has two such chips (identified as host controllers 2442 and 2444) and, as a result, it has two built-in “USB Root Hubs”.

Typically, two of these USB ports are connected to one built-in hub, and the other two are connected to the second built-in hub. Most desktop systems have two USB ports, physically mounted on front side of computer, and two others on the back. Some vendors, such as DELL, place the USB ports connected to the same built-in hub on same side of your computer. Other vendors, like IBM, place the USB-ports connected to different USB-hubs on the same side of the PC. By trying different USB-port combinations for a cascaded NuAmps system, you can find the one in which both NuAmps are connected to same built-in hub. In that case, ACQUIRE will not show the warning, and it will allow you to record EEG.

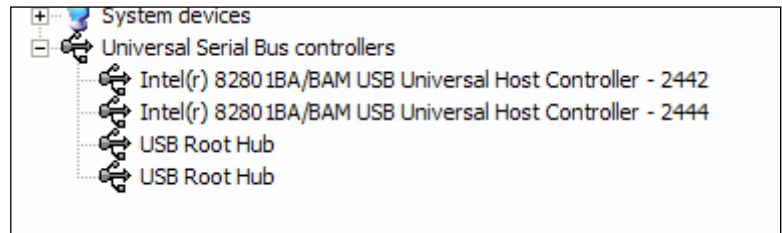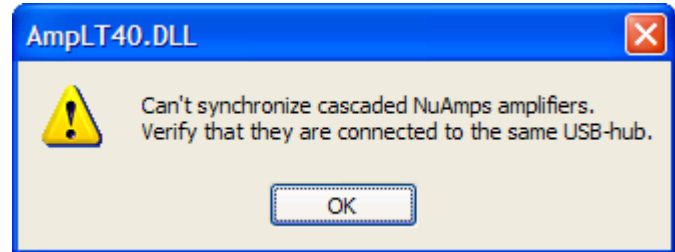

The simplest case is if your computer has two USB ports and only one built-in hub. This is the situation found on some laptop models. In that case, you can plug both NuAmps directly to your computer without the external hub and medical grade power supply.

6. I have SCAN system with single NuAmps unit? Can I upgrade it to cascaded NuAmps system?

Yes, however, this will require the USB-hub and medical grade power supply for it, and it may require a more powerful computer and graphics card. Contact sales@neuro.com for more information.

Troubleshooting

The most common problems while recording from single and cascaded NuAmps amplifier are described below. If your problem is covered in this section, follow the instructions and try to resolve the problem by yourself before calling tech support.

1. NuAmps amplifier is not detected or is not sending signals.

The amplifier may not be detected, or may appear to not be functioning, for several reasons. The amplifier may not be connected to your computer by the USB-cable. The USB-cable may be broken. The Driver may not be installed properly - that can happen if you attached the NuAmps amplifier to your computer before you installed SCAN software. As a result, Windows failed to find the correct driver. Least likely is that there is a hardware failure in the amplifier.

Make sure the amplifier is connected to the computer by the USB-cable. Check the power lid on the amplifier. If the lid is on, proceed to the next step. Otherwise reconnect the USB-cable and check the lid again. If the lid is still off, replace the USB-cable, and try again. If the lid is still off, the amplifier is possibly broken and you need to contact tech support.

If the power lid is on, open the “Device Manager” and verify that you can see the “SCAN LT40” device listed under “Neuroscan USB devices”. If the device was recognized, proceed to the next step. If it is not, reconnect the NuAmps to another USB-port on your computer. If your system does not react to the NuAmps connection, either the NuAmps or the USB-port on your computer is broken. Connect the NuAmps to another PC; if another PC reacts on that connection, then the USB-port on your computer is broken. Otherwise, it is a NuAmps failure and you need to contact tech support.

If a question mark is shown next to the NuAmps icon, or if the NuAmps appears in the “Other devices” section, you need to re-install the driver. Right click on the icon, select “Update driver”, and follow the instructions from the “Installing NuAmps” section above (4. Install the NuAmps Hardware).

For cascaded systems, make sure that the driver is installed properly for *all* attached amplifiers. No question mark for any of “SCAN LT40” devices should be shown in the “Device Manager”. If the device manager shows “SCAN LT40” in the “Neuroscan USB devices” section without a question mark, then the device is recognized by Windows and the driver is installed properly.

Start ACQUIRE and select the “NuAmps 40.AST” file from the “Setup Files\NuAmp” subdirectory. Select “Virtual Ground” as a reference. Attach PocketTrace to the amplifier, turn it on, and set it to generate a sine wave. Start acquisition. The program should start acquiring data and display similar waveforms across all channels. If that does not happen, contact tech support.

2. Acquisition starts but all signals are flat.

There are several possible reasons for flat line signals. The wrong vertical scale may have been selected. Make sure that the vertical scale is selected around 100-200 μV , otherwise you simply will not see a signal. It is possible that the Reference or Ground electrode has not been applied. Make sure you prepare the skin for the Reference and Ground electrodes with maximum care to reduce compound impedance. Verify the physical connection of the Reference and Ground electrodes between the subject's head and amplifier. Use a multi-meter to check the

continuity of the electrode wires. Also make sure that you are using the same Reference that is specified in the software.

The wrong type of electrodes can result in signals that are saturated beyond the working range of the amplifiers. If a high pass filter is applied to a saturated signal, the signal will appear as a flat line. Check the saturation level at the left side of ACQUIRE screen. If the signal is saturated because of depolarization potentials, consider switching to Ag/AgCl electrodes which normally show DC shifts below 10%.

Verify that you have connected the QuikCap to the amplifier. In the case of cascaded systems, make sure the proper adapter is connected to each amplifier.

3. No triggers are detected.

There are several ways in which triggers can be lost. Generally these are easy to resolve, once you determine what the cause is. The simpler causes are with the connections. Make sure the NuAmps-to-STIM cable and stimulation cable are attached appropriately. For a cascaded system, the NuAmps-to-STIM cable should be attached to first amplifier. Verify that all cables are connected tightly and try again. If no triggers are seen, proceed to the next step.

Verify that the stimulation system is powered from the same wall outlet as the SCAN system. Usually this is not a problem, but in some facilities there can be differences in voltages. As a result, the TTL logic is not defined between the two systems, and the triggers are not detected.

If that does not solve the problem, there could be a defective cable (shorted or broken wires). Disconnect the STIM system from the NuAmps-to-STIM cable and test if triggers are working by “NuAmps Setup” dialog, using the “Triggers” page (see Appendix A for more details). Press the “Test” button and short pin 25 (GND) to pins 1-8 and 21-24. As you short the pins, the state should change. If not, then the NuAmps-to-STIM cable is broken, and you need to contact tech support. If the bits themselves are functioning, there may be a problem with the stimulus presentation system (such as a "stuck bit"). Please refer to Appendix A for more detailed testing directions. It is

likely, however, that if the problem is not cable related, there is a hardware failure, which may necessitate the return of some components.


4. The amplifier is not responding.

In most cases, this indicates that the correct driver is not installed. See the Driver Verification section above for more details. If that looks correct, contact techsup@neuro.com for assistance. There are no accessible parts in the NuAmps that can be repaired on-site. A genuine hardware failure will require a return for repair, as will Firmware updates in many cases.

5. Poor USB-performance under XP.

With multiple amplifiers running under XP, there are ways to improve the performance. See the "Tuning Performance of Windows XP Systems" section above for suggestions.

6. Saturation of signals.

Saturation and clipping of signals should be a very rare occurrence because of the broad dynamic input range of the NuAmps amplifiers. Try enabling the  32-bit acquisition option (Overall Parameters \ Amplifiers). This will give the broadest dynamic range. It is better, however, to remove the cause of the drifting. When saturation occurs, it is usually because you are mixing metals among electrodes, especially for the ground and/or reference electrodes (if all channels are affected). It can also occur if the ground/reference electrodes are not well applied (high impedance). See the "DC Electrodes Considerations" section above for more information. It can also occur if the ground/reference electrodes are in need of cleaning or replacement.

7. Missing triggers.

You are seeing some or most of the triggers, but occasionally some are missed. If you are using a stimulus presentation other than STIM, be sure that the duration of the TTL pulses is at least 100 microseconds. Anything shorter than that may be discarded as noise. The optimal TTL pulse width for NuAmps is 500 microseconds.

Assuming the trigger pulses are being sent reliably from your stimulus presentation system, missing triggers can occur if the SCAN and STIM systems are plugged into different power outlets, in which case the TTL logic between systems may be undefined. (In an earlier version of the software, there was a problem in the DLL file that could result in missing triggers; please contact techsup@neuro.com if other potential causes have been eliminated).


8. Increased line noise on some or all channels.


If you see increased 50Hz or 60Hz line noise on all channels, the most likely causes are high impedance (poor application) at the ground or reference electrodes, poor earth ground in the building, or the presence of emanated noise from some nearby equipment (see also the "Notch Filters and Noise" section above). The most likely cause of increased noise on a few channels is high impedance for those channels. If the affected channels are all from the same area (such as the front of the head, decreasing posteriorly), look for a radiating source from some other device in front of the subject (such as an older PC monitor).

The quality of any electrode is related to 1) how the skin was prepared by technician; 2) surface of electrode (Ag/AgCl electrode surface can be damaged); 3) how electro-conductive paste or gel was applied, or if gels with differing conductive properties were mixed. The most common problem is the first one, which is directly related to technician experience. Environmental noise can arise from: 1) proximity of power lines and impulse generating equipment; 2) quality of ground line; 3) crossing of electrode leads by a power line; 4) proximity of the subject's head to noise emitting equipment, e.g. computer screen. In addition, if you place an electrode with high impedance ($> 15\text{K}\Omega$), NuAmps will still record EEG because of its high input resistors ($=100\text{M}\Omega$), but then any patient movement or movement around him will generate an "antenna" effect.

Appendix A: Trigger Port Interface

The method of triggering between the SCAN and STIM systems is to connect the Parallel I/O ports on the back of the SCAN and STIM interfaces by means of a STIM-to-SCAN cable. This cable is provided to customers who purchase STIM systems. A STIM-to-SCAN cable is a grey, 25-line cable with a female connector at the end marked SCAN and a male connector at the end marked STIM. A parallel-to-serial cable is provided with NuAmps to allow input into the NuAmps unit. Below are diagrams of the STIM-to-SCAN cable connections.

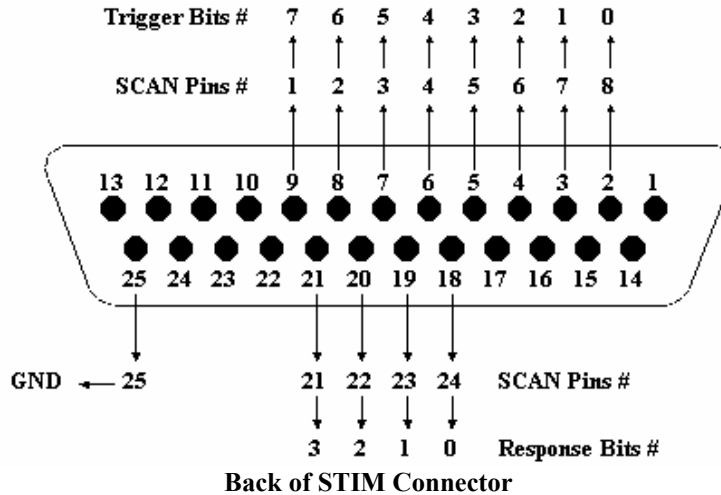

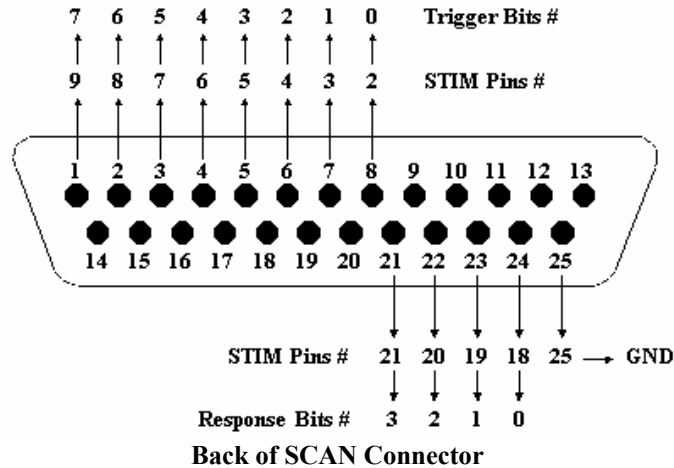

Response pad pins

On the back of the STIM connector the response pad lines are 21, 20, 19, 18, carrying response bits 3, 2, 1 and 0. On the back of the SCAN connector the response pad lines are 21, 22, 23, and 24, carrying response bits 3, 2, 1 and 0.

Matching Port Logic

The logic used to trigger your system can be either positive or negative. Positive logic is defined as a transition from the zero state (ground) to a one state (5 volts TTL). Negative logic is defined as a transition from the one state to the zero state. The NuAmps employs positive logic on the stimulus port and negative logic on the stimpad (response port). *Note - the stimpad lines must be held to positive levels during the resting state or the device will not respond. If not connected these lines are 'pulled-up' to high levels.*

Because NuAmps uses positive logic on the stimulus port the numbers received by the trigger port from older 100kHz based STIM systems (identified by the two large ribbon cables exiting the back of the computer labeled J6 and J11) need to be inverted back into positive logic to match the NuAmps. This is accomplished by the **Use inverted values** option in Acquire, under Setup / Edit / Triggers (see Acquire manual). In general, if trigger codes do not match (i.e., a 1 generates a 254) then the **Use inverted values** option should be used.

Troubleshooting

If you encounter problems with triggering, there are a number of procedures that you can follow to determine the cause of the problem. Most triggering problems are simple to correct. These are described completely in an Appendix at the end of the Installation/Orientation manual. The sources of the problem can be on the STIM side, the NuAmps side, or the connecting STIM-to-SCAN cable. Only those problems relevant to the SCAN side with NuAmps will be discussed here (refer to Appendix mentioned above for other potential causes).

The NuAmps Setup option, under Edit, has the Trigger Test tab. This reads the status of the bits at the trigger port on the NuAmps. It provides a quick means to see what the problem is and where the source is.

The Decimal column will show the type code numbers that are sent from STIM. The Binary column will show the bits as they change. Bit 0 is on the far right side of the display, and bit 7 is on the far left side.

Stimulus. The Stimulus fields display the stimulus trigger type codes received. With STIM running, you will see brief flashes in the Stimulus decimal and binary fields with each trigger received. If you increase the pulse duration in the STIM software to, for example, .1s, you will see the triggers somewhat more

NuAmps-to-STIM cable testing

	Decimal	Binary	
Stimulus	0	00000000	Test
Last Stim	0	00000000	
Response	0	00000000	
Counter	0		
Hold Value	0		

clearly. The type code number will flash on the decimal field, and the activated bits will flash on the binary field.

Last Stim. The Last Stim fields will show the type code in decimal and binary of the most recent trigger.

Response. The Response decimal field will show the inverted value of the response pad triggers (i.e., 255). Pressing the response pad buttons should show decimal values of 254, 253, 251, and 247, corresponding to buttons 1-4, and type codes of 1, 2, 4 and 8, respectively. The binary field will be all 1's until a button is pressed, then the corresponding 1 will become a 0 (uses the 4 columns on the right side of the display).

Counter. The Counter will show the accumulated number of stimulus and response triggers received.

Hold Value. The Hold (bit pattern) value determines what value of the signal at the trigger port will initiate a sweep; a value other than the Hold value will trigger acquisition. With new STIM systems, the Hold value should be zero. If you see a different value, please contact Technical Support. If you have an older LabMaster STIM system, and do not see any triggers with a Hold value of zero, you need to invert the trigger codes. Do this by entering 255 for the Hold value from the keyboard (you will need to click the *Use inverted values* field under Edit/Overall Parameters/Triggers). If you still do not see triggers, contact Technical Support.

With the STIM system turned on (but not sending triggers) and the stim-to-scan cable connected, go to Trigger Test, and click its Test button. If everything is functioning normally, you should see:

In its “resting” state, note that the Stimulus bits are at zero, and the four Response bits are held high (all 1’s). Summing the 4 response bits gives the 15 Decimal value. The Hold Value should be 0.

If you see, for example, one (or more) of the stimulus bits held high (1), that will cause all other incoming triggers to be ignored (no triggers). A workaround is to take whatever number it shows in the Stimulus field, and enter that as the Hold value*. For example, if the 4th bit is stuck (8), enter 8 as the Hold value. (The events seen in ACQUIRE will be altered accordingly - this is a temporary workaround until the real problem is repaired). The cause could be a problem in the STIM box, the stim-to-scan cable, or the trigger input circuitry in the NuAmps.

	Decimal	Binary
Stimulus	<input type="text" value="0"/>	<input type="text" value="00000000"/>
Last Stim	<input type="text" value="0"/>	<input type="text" value="00000000"/>
Response	<input type="text" value="15"/>	<input type="text" value="00001111"/>
Counter	<input type="text" value="1"/>	
Hold Value	<input type="text" value="0"/>	

With other stimulus systems connected, you might see all response bits at 0. They all need to be high in the resting state. Then it is a question of whether you want to record responses or not. If you do not want responses, then you should not plug anything into those pins on the SCAN side of the trigger connector (pins 17-24). The natural resting state of the NuAmps is high, so they will be OK. If you do want responses, they must use inverted logic, where the resting state is high, and the trigger pulse goes to zero.

Now start a program in STIM that sends stimulus triggers, and/or press the buttons on the response pad. You should see the corresponding bits (the 1's and 0's), toggle on and off. These will be very brief - increase the duration of the trigger pulse in STIM to make the toggling easier to see.

To help isolate the cause of abnormal bits, disconnect the STIM-to-SCAN cable from the back of NuAmps (and Start Test trigger port). You should see:

	Decimal	Binary
Stimulus	255	11111111
Last Stim	255	11111111
Response	255	11111111
Counter	1	
Hold Value	0	

All of the bits should be high. If any are at zero, that points to a problem in the NuAmps (and a probable return to El Paso). If this looks normal, connect the STIM-to-SCAN cable to the NuAmps, and disconnect it from the stim box. It should still look like the picture above. If it doesn't, there is a problem in the cable. If it still looks normal, then the problem is likely on the STIM side.

When you are testing for triggers in ACQUIRE, using the Single Window display to see the triggers in a continuous file, be sure you have NOT enabled the "Use inverted values" field in the setup file (look under Edit, Overall Parameters, Triggers).

Appendix B: NuAmps Electrode Starter Kit

Your system will arrive with an Electrode Starter kit, and/or a NuAmps Quik-Cap(s)[™], if you purchased them. Application of the electrodes will vary depending on which you have.

The NuAmps Starter Electrode Kit contains 40 electrodes and a container of Ten20 paste, with some additional accessories. These NuAmps electrodes feature:

- Silver / silver chloride electrodes are standard, but, by special order, you may have received sintered silver or gold electrodes.
- All discs are 10 mm (0.4 in) in diameter and come with 1.5 m color-coded wires. That length ensures a comfortable recording environment.
- The electrodes have touch proof protectors, which prevent the wire from breaking down while bending in any direction. This dramatically extends the life of the electrodes.
- EEG wires are coated with Teflon to resist tangling. Electrodes that are jumbled together easily pull free.
- Electrodes are cupped, with the hole at the peak, to facilitate better contact with electrode paste. The hole is used to refill conductant during prolonged recordings.




Ear clip electrodes are not necessary, and any electrode can be placed at the reference point.


We are assuming that registered EEG technicians will be recording the EEGs, and the information below should augment your existing knowledge.

With individual electrodes, you can place them anywhere, although the International "10-20" system is used for labeling in NuAmps. We recommend you adhere to the 10-20 (or modified) system to avoid confusion. The caps follow the modified 10-20 system, with its additional electrodes.

The basic 10-20 system is as follows. Electrodes are placed over the various lobes of the brain, and the electrode labels identify the lobes over which they are positioned. Odd numbers are from the left side of the head, even numbers are from the right side, and electrodes down the middle have a "z". (While there is no central lobe, the label C is used for convenience). For example, the F3 electrode is from the frontal lobe on the left side. Cz, or Czcd, is placed at the vertex at the top of the head. Precise placements use percentages of distances measured from landmarks on the head - the nasion (the dented junction of the nose and the forehead), and the inion (the bony protuberance at the back of the head).


Diagrammatic placement of the standard 10-20 system electrodes is as follows.

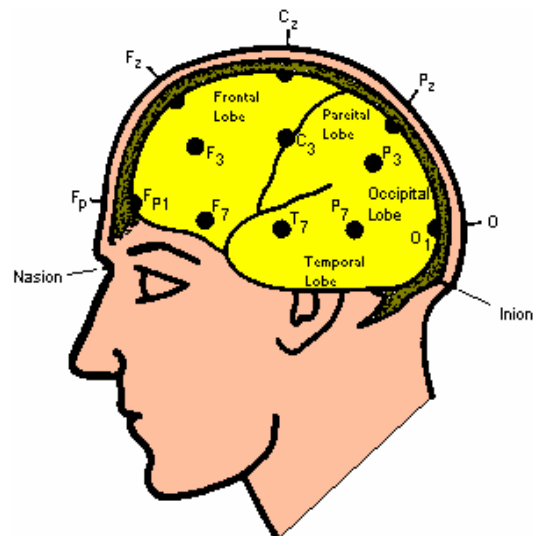

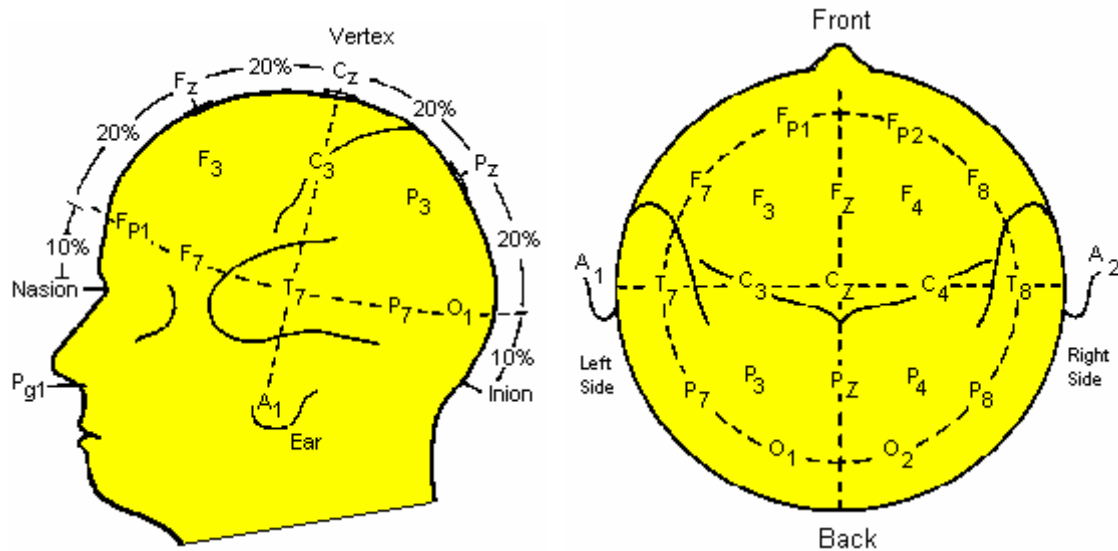

The modified 10-20 system uses these same electrode positions, as well as additional ones that fall in between these positions. For example, FC3 lies halfway between F3 and C3. An approximation of the modified system is below.

Appropriate spatial sampling of potential fields requires regular placement of electrodes at sufficiently close intervals across the scalp. The figure above shows 19 labeled electrode positions having 20% spacing along standard lines of measurement, plus additional electrodes, FPZ and OZ, having 10% separation from nearest neighbors along a coronal line (FP1 and FP2 are the nearest 10% neighbors to FPZ, and O1 and O2 are nearest to OZ). The 20% interval separates nearest neighboring electrodes by 5-7 cm in adults. Although this is usually adequate for routine EEG, more precise localization, as needed in topographic mapping or delineation of epileptogenic foci, requires placing intermediate electrodes to reduce spacing to 10%, or about 3 cm.

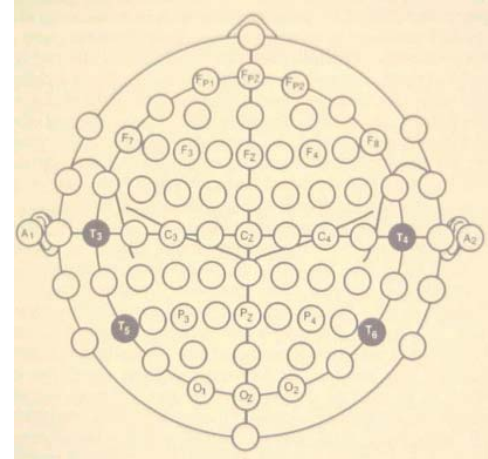

Electrode Application. Because NuAmps amplifiers are true DC amplifiers, you will find there are some differences in the application of electrodes.

- If you are using gold electrodes, they must be presoaked in water for about 2 hours before use.
- Do NOT mix metals across electrodes, e.g., do not use silver electrodes for the active channels and gold electrodes for the reference(s).
- Do NOT use one type of conductant at some sites, and a different type at other sites.
- Do NOT vary the skin-electrode interface, that is, if some of the electrodes do not actually touch the scalp/skin, as is the case with electrode caps having recessed electrodes, then all electrodes must be off of the scalp/skin. In the case with caps, this is done by placing a double-sided circle or V-shaped strips of tape

between the electrode and skin for the reference(s) and EOG leads. The tape will also help hold the electrodes in place.

Failure to follow these directions can create DC battery potentials, which may, in turn, lead to clipping of the signals or complete saturation of the channels. The EEG will be either distorted or lost completely.

When using individual electrodes. Prior to placing an electrode, the scalp must be prepared by rubbing it vigorously with alcohol or a skin preparing agent. This removes dirt and oil from the electrode site and lowers the impedance. Obviously, excessive cleaning can be irritating to the patient. Similarly, some patients are sensitive to application of paste containing salt solutions or bentonite. Several varieties of electrode attachment media are available, including sodium chloride pastes or gels, conducting sponges, and other specialized electrolytes. Electrodes may be held in place by viscous gels, by mechanical restrictions (tapes, V-shape strips, bands or rubber caps), or by collodion.

If you use an abrasive paste, it is recommended that you use it on all sites. With DC amplifiers, it is important to treat all of the electrode contacts the same. The same type of electrodes (metal), the same conductant, and the same interface should be used with all electrodes. Mixing metals, or conductants, or skin/scalp interfaces can create DC battery potentials that may then lead to clipping or saturation and loss of EEG signals.

For prolonged recording, collodion is a glue formulated from pyroxylin in ether, alcohol, or camphor; it is liquid when applied, but dries to a strong adhesive within minutes. It is suitable for patients who cannot keep still, or who will need electrodes in place for more than a few hours. Pyroxylin is an element of gunpowder, and ether is highly flammable. Neither flame, nor excessive heat, nor pure oxygen should be in the vicinity of collodion applications. One must remember that collodion is not an adequate conductive medium (do not let it get between the electrode and the scalp). Conductive gel must be injected into the cup electrodes and refreshed periodically. A blunt needle is usually employed for this task. Collodion is removed with acetone scrubs.

With NuAmps, you MUST place the GROUND electrode first, then place the REFERENCE electrode(s) in order to see valid impedances as other electrodes are attached. Pay particular attention to attaching the ground and reference electrodes - attach them as securely as possible. Then place the active EEG electrodes.

Typically, after placing the GROUND and REFERENCE(S), you would place the CZ electrode next. Measure the distance between the Inion and Nasion using a measuring tape (in cm's). Calculate half that distance. Then measure from just in front of one ear to the same place in front of the other ear (positioning the tape over the top of the head). Ask the patient to open and close his mouth - the "dent" you feel in front of the ear is the point to take the measurements from. Calculate half of that distance. The intersection between those two imaginary lines - the half way point of each - is the location of CZ. The other placements are based on percentages (usually 10% or 20%, hence the 10-20 system) of those measurements. For example, O1 is 10% up from the Inion, and 10% over to the left. (If you use a wax pencil to mark the locations, be sure to remove the wax when you scrub the area. Wax is not a good conductor and the impedances could be affected).

P3 is 30% up from the Inion, and 20% over to the left. T7 is 10% up from the left ear measuring point (on a line through CZ). C3 is halfway between T7 and CZ. P7 is halfway between T7 and O1, and so forth.

The X1 electrode (for example) is placed above the left eye (VEOG1+), directly above the pupil (looking forward), and between the eye brow and Fp1. X2 (VEOG1-) is placed on the skin directly below the pupil, on the cheek bone, about 2cm below the eye. It is typically not necessary to place VEOG electrodes for both eyes, since the activity from the eyes is highly correlated. X3 (HEOG1) may be placed about 2cm to the left of the left eye outer canthus (corner of the eye). X4 (HEOGr) is placed about 2 cm to the right of the right eye outer canthus.

After the recording, clean the electrodes gently, taking care not to scratch or pit the electrode surface. This can also lead to DC problems. Disinfect the electrodes according to safety protocols established in your facility. For example, soak the electrodes in a disinfecting solution (such as Abcocide, Metricide, or Betadine) in a concentration of 1 part

disinfectant to 4 parts water, for 15-30 minutes, then dry them before the next use. If full strength disinfectant is desired, the electrodes should not be left in solution for more than 10 minutes.

For maximum electrode life, do not expose reusable surface electrodes to prolonged soaking or to cleaning with corrosive agents or bleach.

NOTE: electrodes should be cleaned and disinfected, but not sterilized because they are used externally.

Neuroscan recommends that you inspect electrodes for extensive wear or damage prior the use. If the electrode wire, termination, or surface area is worn so as to impair performance, the electrode should be discarded. Neuroscan is not responsible for injury, infection, or other damage resulting from improper electrode preparation or use.

When using electrode caps. To apply the reference electrodes to A1 and A2, for example, first clean the ear lobes with an alcohol prep pad, and let the alcohol evaporate. If you then use a preparation paste, gently scrub the area with that. In this case, you should use the same paste at all of the electrode sites. With electrode caps, this can be accomplished by using the wooden end of a Q-tip that is dipped in the paste, inserted through the electrode hole, and then gently rubbed on the scalp. The purpose is to make sure that the chemicals that are present between the electrodes and the skin/scalp are the same at all channels (to avoid potential DC battery potentials). If you find that your particular paste does not affect DC levels, you can safely ignore this precaution.

Place either a double-sided circle or V-shaped strip of tape on the earlobe or skin surface, and then place the reference electrode(s) on the tape. Inject the conductant, as described below, using a gently abrasive movement with the blunt syringe tip on the earlobe. Adding a piece of surgical tape over the top of the electrode and around the ear lobe will also help hold it in place. (If you need to add more gel, you can carefully poke the syringe through the tape). Attach any other loose EOG electrodes in the same fashion.

Quik-Cap Application. Your Quik-Cap has been designed with a highly elastic fabric that provides a uniform fit over a wide range of head size and shape variability. Quik-Caps currently use silver/silver chloride electrodes (gold and sintered silver are other options). The electrodes are carried in soft rubber holders for added wearer comfort.

When donning the cap, locate it loosely on the head after positioning the earlobe or mastoid electrodes first, using recommended skin preparation and electrolyte application techniques. Next pull the cap down firmly over the head, attach the chinstrap and ensure that the FP1 and FP2 electrodes are approximately 4cm above an imaginary line extended out from the Nasion. Also, check that the rear of the cap is stretched down the neck while holding the front of the cap in place, and then tighten the chinstrap slightly to fix the cap. Finally, place the VEOG monitoring electrodes, as described above for individual electrode placements.

Now that the cap is positioned, attach the wiring with tape or other fastening means (such as a Velcro® strap) onto the shoulder of the wearer. Ensure that there is adequate slack in the cable to allow head movement without pulling the cable. The final step is to insert a blunted syringe through the hole in each electrode down to the scalp surface. Lightly rock the syringe while moving it in a circular fashion and inject a small amount of EEG electrolyte. Make sure that you lift the syringe while injecting the electrolyte to guarantee that there is an unbroken path of electrolyte through the hair to the electrode disk (with no air bubbles). Be sure not to inject too much electrolyte because this could cause neighboring electrodes to be shorted together. Fill all electrodes in the cap before measuring impedances at any site except for the reference location(s) to allow the impedance values to stabilize. In order to minimize the amount of skin abrasion and increase the comfort of the wearer, we recommend using Quik-Gel EEG electrolyte with all cap systems. (NOTE: DO NOT COMBINE OTHER TYPES OF ELECTRODE MATERIALS OR USE MORE THAN ONE TYPE OF ELECTROLYTE WHEN USING THE QUIK-CAP.)

Using Quik-Gel. Quik-Gel is unique in its formulation and requires a few simple techniques to reduce your EEG electrode application time. Quik-Gel can absorb moisture from the scalp's dead skin barrier to rapidly lower electrode impedance. Therefore, little or no skin abrading is necessary to attain the same result as with other brands of EEG electrolyte.

Follow accepted selection and preparation techniques for reference, ground and EOG electrode locations. Inject a small amount of Quik-Gel into each electrode cup and apply to the skin using normal methods, being careful not to bridge the electrolyte between neighboring electrodes. For application into thick hair, press firmly on the electrode housing, insert the syringe tip and gently abrade as mentioned above. Note the following tips.

- Ensure that the syringe tip is lifted through the hair while Quik-Gel is being injected in order to ensure an unbroken contact between the scalp and electrode surface (with no air pockets).
- Apply electrolyte to all electrodes before measuring electrode impedance. In this manner, Quik-Gel will work to reduce the input impedance to acceptable levels necessitating less skin abrasion and patient discomfort.
- Use an elastic gauze size 4, 5 or 6 to cover the cap. You will find a more uniform pressure of the electrodes and a more stable impedance reading.
- When exposed to air, Quik-Gel will begin to crystallize. Place a glass of water next to your installation and when the syringe is not in use place the tip in the glass.
- KEEP CONTAINER TIGHTLY CLOSED WHEN NOT IN USE
- Do NOT refrigerate Quik-Gel.
- Distilled water can be used to thin the mixture. Be sure to shake well to ensure uniform consistency.

Impedance testing. The steps to follow have been described above, and only the main points to remember are repeated here. When using individual electrodes, you *MUST securely* place the GROUND electrode first, and then the REFERENCE(S) in order to see valid impedances. Pay particular attention to the ground and reference electrodes.

If you are using Quik-Gel, it is a good idea to place all of the electrodes first before looking at the impedances. Quik-Gel minimizes the amount of abrading that is necessary when you allow it to soak in for a few minutes.

Impedances may be read directly from the Impedance display, and they should be equal to or less than 10 kOhms. They can also be read from the display on the face of the amplifiers. The ground impedance is not displayed, but the reference impedance is, and you should get the reference impedance(s) as low as possible. If they are high, all of the channels will be high. NuAmps will record with impedances up to 200 kOhms (although the higher the impedance, the greater the likelihood of noise on the recording).

Quik-Cap Cleaning and Disinfecting Instructions. The following method is recommended for cleaning and disinfecting your Quik-Cap. After each use, remove all tape, markers, and electrode pads from the cap and unhook any wire harness shoulder restraint. Unplug the wire harness from the amplifier equipment and remove the cap from the wearer. Take the cap to a sink and rinse off any electrolyte on the outside of the cap and electrodes. Next, turn the cap inside out to expose the electrodes that are down inside the rubber holders. Run a tight stream of very warm or hot water from the faucet and use it like a drill to remove the electrolyte from the holder cavity. This should remove the electrolyte completely if you are using Quik-Gel, if not, and the electrolyte has oils in it (like ECI-Gel), you will need to use cotton swabs to help remove the electrolyte residue from the electrode surface while the water is running. Do this carefully to avoid scratching or pitting the surface of the electrodes. Next, turn the cap back out and use the water to remove any additional electrolyte from the outside of the cap.

When you are satisfied that nearly all electrolyte has been removed, place the cap and wire harness into a warm water and soap bath being careful not to get the connector at the end of the harness wet. The bath should contain 4 quarts of water with 1 or 2 ounces of mild dish soap (Ivory™, Dove™, etc.). Let the Quik-Cap sit in the bath for approximately 1/2 hour. Remove and rinse the Quik-Cap with warm water and place it into a disinfecting solution (such as Abcocide, Metricide, or Betadine) in a concentration of 1 part disinfectant to 4 parts water. If this disinfectant method is selected, leave the cap in the solution for 15 to 30 minutes. If full strength disinfectant is desired, the cap should not be left in

solution for more than 10 minutes. After disinfectant, rinse the cap thoroughly and hang it up to dry. For faster drying, the cap can be set in front of a fan, however, direct heat sources should not be used to dry the cap.

(Note: compressed air can be used to blow moisture from the electrodes and wire entry points and is very useful for extending the life of your Quik-Cap).

Appendix C: PocketTrace-2 Signal Generator

The PocketTrace and PocketTrace-2 (PT-2) units are small and extremely useful devices for generating signals that may be sent to the NuAmps or other EEG/ECG amplifiers. The signals can be used for calibrating the system, and as a means to test the basic functionality of the system without having a live subject connected. Together with the PT-2 Adapter, which has attenuators from 2000:1 to 20000:1, the device will produce a variety of signals with varying frequencies and amplitudes (within the physiologically normal range for EEG and ECG). A TTL pulse is generated with each repeated signal to test trigger functions.

PT-2 is a modified version of the original PT-1 model. PT-2 has 4 control buttons (rather than three in the PT-1 model), as shown on the picture below. It provides the same functionality as PT-1, with better precision ($\pm 1.2\%$ instead of $\pm 3.5\%$), but is controlled differently and has additional information on the LCD. If you need information for the PT-1 device, please contact techsup@neuro.com. The information below is for PT-2.



Note: When using PT-2 in ACQUIRE, you must set the Reference to Virtual Ground. (From the NuAmps Setup Editor, click the Edit button, then select the Reference tab, and then select the Virtual Ground option). Return the Reference option to your original choice when you are finished with PT-2.

Cable Connections

Please refer to the diagram on the next page showing PocketTrace connected to the NuAmps amplifiers. In this configuration, the signal from PocketTrace will be sent to all channels in NuAmps.


Connecting the cable(s) to NuAmps. The supplied cables have a large connector on one end (DSUB37-M connector), and 7 cables coming out from the back of the large connector. Six of the wires end in touch-proof connectors. One is labeled GND (ground), one is labeled ECG (electrocardiogram), and the other 4 are not

labeled. The seventh wire is longer and larger, and it ends in a mini-phone connector. You can take the output from PocketTrace and send it through all the channels, or you can selectively choose up to 5 input channels.

Sending the signal to all channels. To send the signal generated by PocketTrace to all channels, simply plug the large connector into the CAP connector on the NuAmps amplifiers. Then connect the longer (and larger) wire into the SIGNAL Output jack on the side of PocketTrace, and connect the smaller, loose wires as shown in the diagram on the next page.

Sending the signals to selected channels. Connect the GND lead to either GND jack on the amplifier unit. The ECG lead differs from the 4 unlabeled channels in that its output has a much larger amplitude. Plug it into any desired channel (or leave it disconnected). Plug the remaining leads into any other desired channels (or leave them disconnected). Connect the longer (and larger) wire into the SIGNAL Output jack on the side of PocketTrace. Do NOT connect the large connector to the CAP connector on the side of the NuAmps amplifiers.

Connecting the TTL pulse cable (optional). With some of the Modes (described below), a TTL pulse is sent out through the SYNCH Output jack. This is used to simulate trigger inputs into NuAmps (or other device). The output is a positive logic (resting at 0Vs), TTL compatible (5V) pulse with a 10ms duration. This feature currently has no use with NuAmps systems (as TTL pulse triggering has not currently been implemented). NuAmps users should ignore this option for the time being.

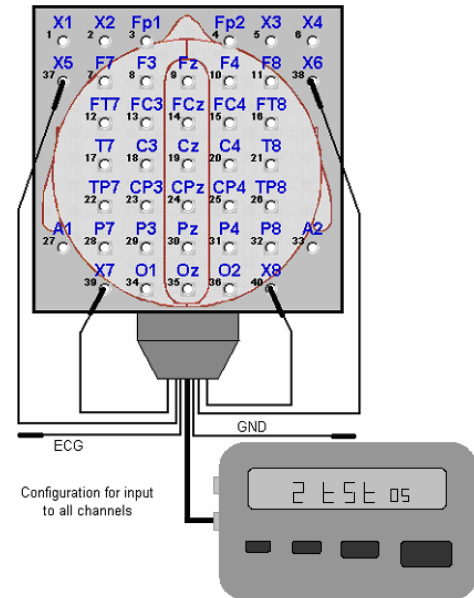

Operations. The follow information explains the basic operations of PocketTrace.


Turning the device on/off. Press the left-most button (with the red dot) to turn on the device. Press the same button again to turn in off. The device will turn off automatically if you do not press any button for 10 minutes. If you press the “+” button and then turn on the device, it will operate until the battery is discharged, or until you turn it off. With this mode of operation, the “globe” symbol is shown in the right upper corner.

Blinking indicator. If the LCD blinks, you need to change battery.

Signal type selection. Use the “Mode” button to change the type of generated signal.

Selecting the frequency of the output signal. The Increase (+) and Decrease (-) buttons are used to change the *frequency* of the output signal. This frequency is shown on the right field of LCD. A single press of the Increase button will increase the frequency by one step. If this button is pressed for several seconds, the frequency will increase rapidly and continuously. The Decrease button operates in same way to decrease the frequency.

Selecting the amplitude of the output signal. To increase the amplitude, press the “+” button and while keeping it pressed, press the “Mode” button. To decrease the amplitude use the “-“ and “Mode” buttons. The

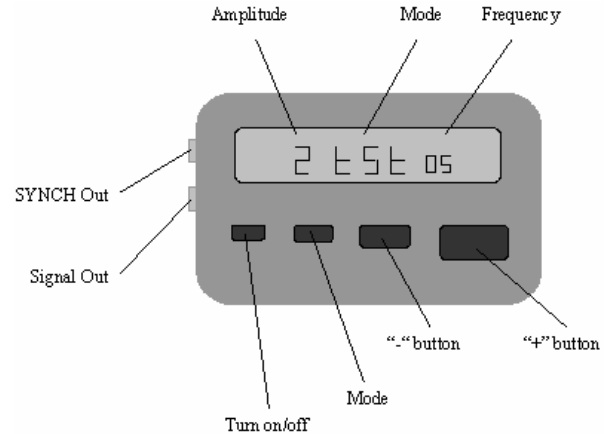

amplitude of the output signal is indicated by the symbol to the left of the signal type. The device provides 3 amplitudes:

Full scale ($\pm 1.25\text{V} \pm 1.2\%$) – no symbol at the left is shown

Half of scale ($\pm 0.65\text{V} \pm 1.2\%$) – symbol at the left is “1”

Quarter of scale ($\pm 0.31\text{V} \pm 1.2\%$) – symbol at the left is “2”


With the half-scale and quarter-scale modes, the signal has DC components. In "Square pulses" mode with a frequency setting of 0Hz (described below), the signal amplitude is 0, and the DC component varies as shown:

Amplitude Indicator	DC shift
None	-1.23V, $\pm 3.5\%$
1	0V, $\pm 45\text{mV}$
2	+1.23V $\pm 3.5\%$


Signal output jacks. The signals from all modes in PocketTrace are sent out from the SIGNAL Output jack on the side of the device. With some modes, TTL pulses (positive logic, voltage pulses) are sent out from the SYNCH Output jack on the side of the unit.

Operational modes. There are 6 different operational modes: a square wave generator, a sine wave generator, 3 visual evoked response simulations, and an ECG simulation.


“Square pulses” mode. In this mode the device produces rectangular square wave pulses in a frequency range of 0.5-99 Hz.

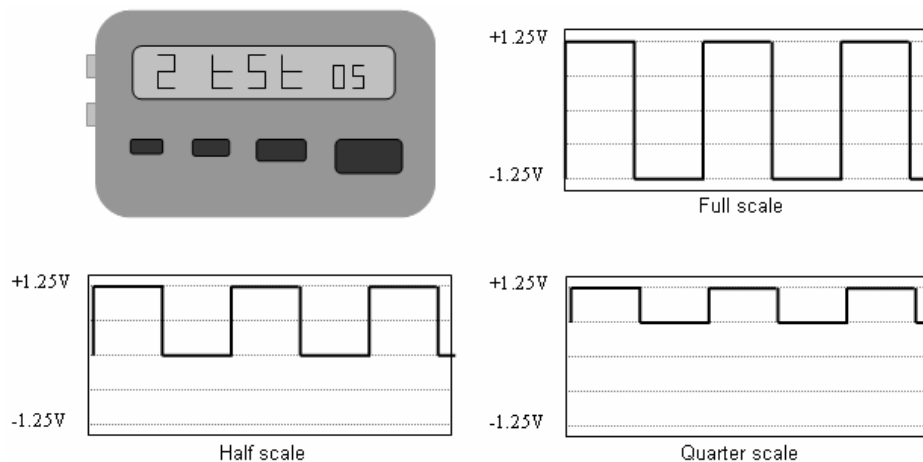


In ACQUIRE, you will see signals similar to the following (High and Low Pass filters OFF):

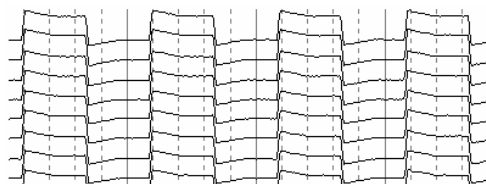

“Sinusoidal signal” mode. This mode is used to generate a sinusoidal signal in the 0.5-90 Hz frequency range.

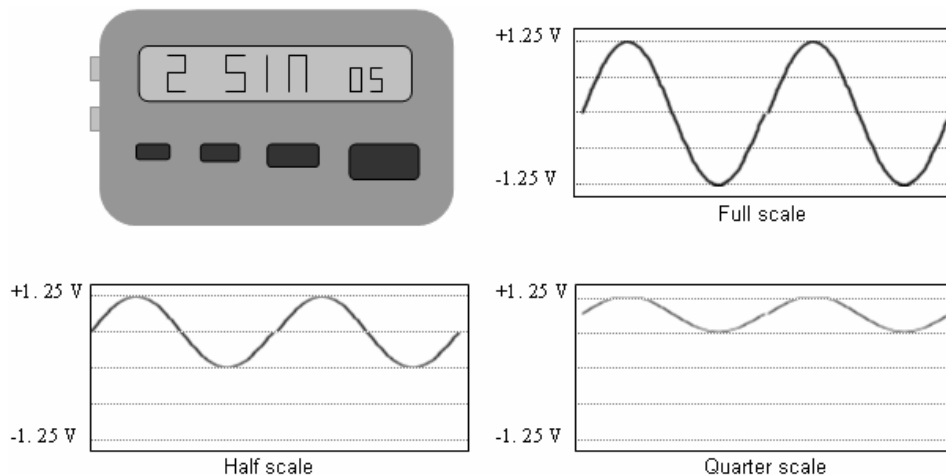


In ACQUIRE, the signals will look similar to the following:

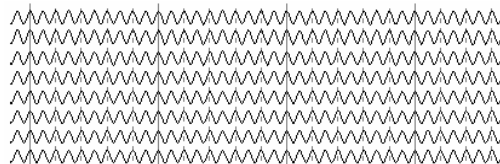


“Special wave-EE1” mode. In this mode the device produces a simulated visual evoked response to a *checkerboard reversal* pattern. The signal consists of a series of impulses, which appear with a frequency range


from 0.5 to 1.6 Hz. A new impulse can be generated even if the previous impulse has not finished (the impulses may overlap). Simultaneously, with every new impulse generation, a TTL trigger pulse (positive logic) is sent from the SYNCH jack. *Note: the frequency of impulses is expressed in Hz times 10 (for example, display “15” corresponds to 1.5 Hz).*

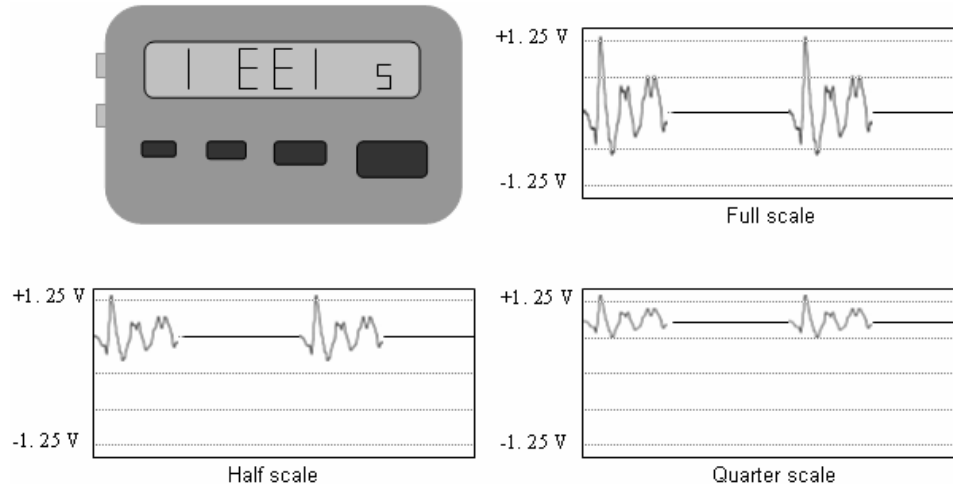


In ACQUIRE, the signals will look similar to the following:

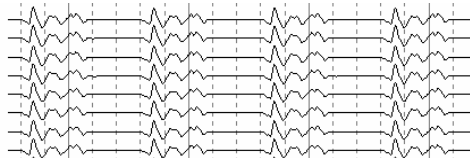





“Special wave-EE2” mode. In this mode the device simulates a visual evoked response to *flash* stimulation. The signal consists of a series of impulses in the frequency range of 0.5-1.6 Hz. A new impulse can be generated even if the previous impulse has not finished (the impulses may overlap). Simultaneously, with every new impulse generation, a TTL trigger pulse (positive logic) is sent from the SYNCH Output jack. *Note: the frequency of impulses is expressed in Hz times 10 (for example, display “15” corresponds to 1.5 Hz).*

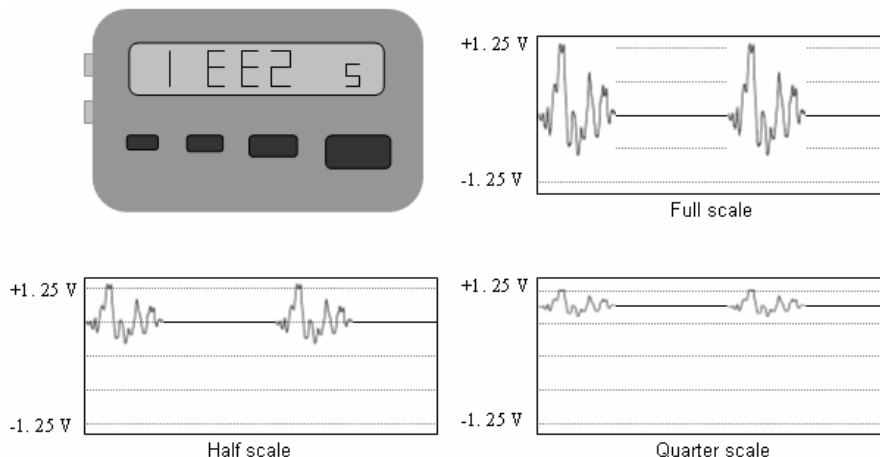


In ACQUIRE the signals will look similar to the following:

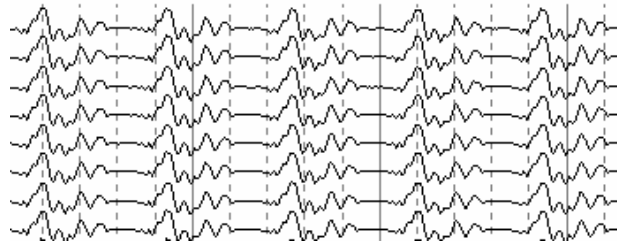


“Special wave-EE3” mode. In this mode the device produces alternating “special wave-1” and “special wave-2” signals from the SIGNAL Output jack. The signals have a frequency 0.5-1.6 Hz. *Note: the frequency of impulses is expressed in Hz times 10 (for example, display “15” corresponds to 1.5 Hz).*

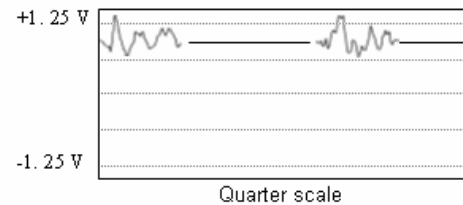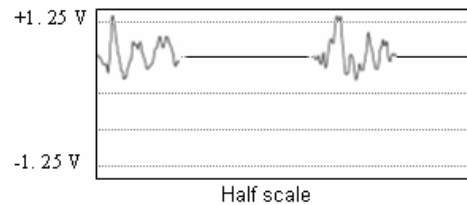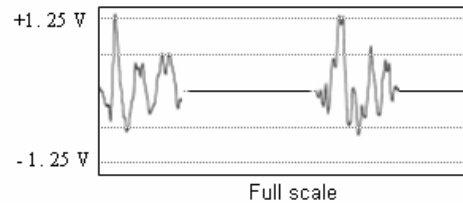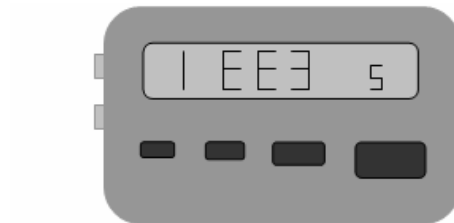

In ACQUIRE, the signals will appear similar to the following display:

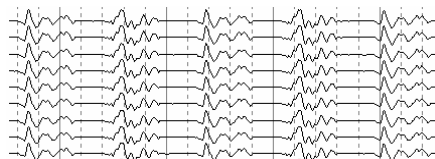

ECG Mode. This mode simulates ECG activity. The waveforms appear within a frequency range of 0.5 to 1.6Hz, with three amplitude levels. Simultaneously, with every new impulse generation, a TTL trigger pulse (positive logic) is sent from the SYNCH output jack.

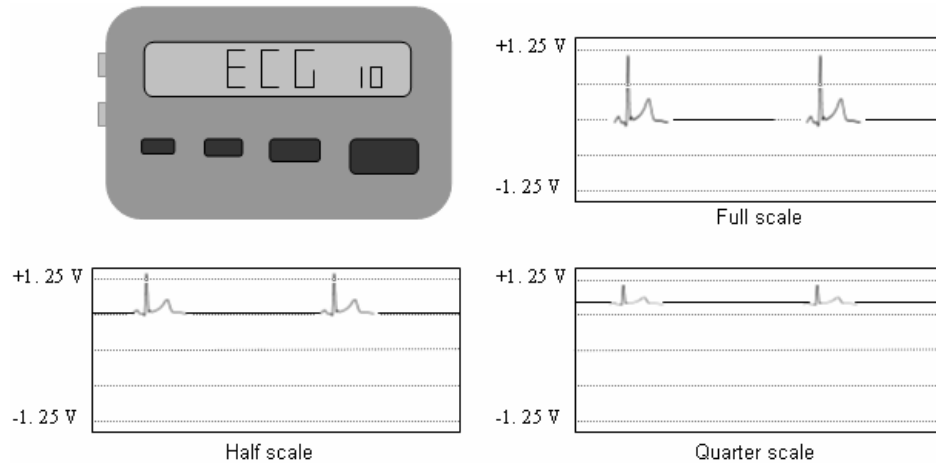


In ACQUIRE, the signals will look similar to the following:



Technical specifications


Description	Value
Peak-to-peak signal voltage on SIGNAL output	± 1.23 (0.62,0.31) V $\pm 3.5\%$
Relative error of signal frequency	Not more than 0.5%
Output impedance for SIGNAL output	Not more than 100 Ω
Synchro-pulse level on SYNCH output	TTL compatible
Duration of synchro-pulse	10 ms $\pm 10\%$
Power supply	One AAA-size battery
Duration of operation for one battery ("Alkaline" type)	Not less than 45 hours

Pin-out of SIGNAL connector

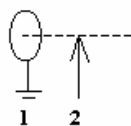

- 1 – common ground
- 2 – output signal

Pin-out of SYNCH connector

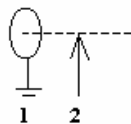


- 1 – common ground
- 2 – TTL trigger pulse

Appendix D: Cap Connector Pinout

The pinout for the cap connector is as shown:

X5 (FT9), X6 (FT10), X7 (PO1), X8 (PO2) electrodes are not wired to DSUB37 but can be accessed through front-panel individual touch-proof connector.

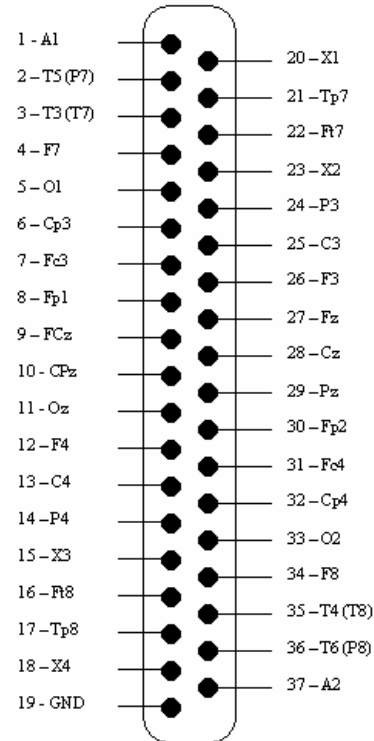


Appendix E: Impedance Measurement

The following section explains how NuAmps conducts impedance measurement.

Scheme of input cascades during impedance measurement

The picture below depicts impedance measurement on a subject's head (outline by dashed circle) with 3 electrodes attached – active (**Ch**), reference (**Ch_ref**) and ground (**GND**) electrodes:

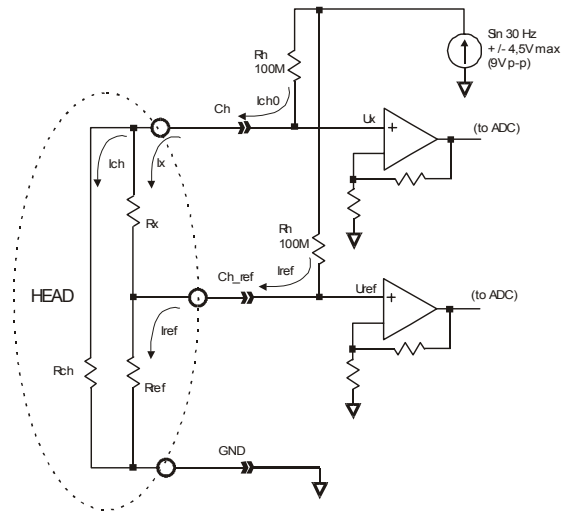

A 30Hz sine signal is injected to the head at all channels, including the reference, through their 100M Ohm input resistors. As a result, each electrode injects a 30Hz harmonic current with 45nA maximum amplitude ($4.5V / 100M\Omega$). Current from all channels flows back to the amplifier through the GND electrode. Therefore, the maximum summary current through GND is equal to the number of active electrodes, plus one reference electrode, multiplied by 45nA.

Denote the current injected through the **Ch** active channel as **Ich0**. This current splits into 2 parts: **Ich0 = Ich + Ix**, where **Ix** is the current through an unknown **Rx**.

Ch can be any channel and **Ch_ref** is the reference channel selected for impedance measuring. Impedance is measured and calculated between **Ch_ref** and the other channels. No special commutation circuitry is used for **Ch_ref**. It is specified by the software and may be any channel. Usage of **Ch_ref** also

minimizes the effect of multiple current sources and better blocks the common-mode voltages (by measuring $U_x - U_{ref}$).

Equivalent scheme for transhead impedance calculation

As impedance inside the head is distributed, finding **R_x** in the case of many active electrodes is a very complicated mathematical problem. Consider the simplified equivalent schematic, where **R_{ref}** = const = 0. As a result, **R_{ch}** transforms to **R_x**:

According to this schematic: $U_{xd}/R_x = U_g/(R_h + R_x)$

After reorganization: $R_x = R_h * U_{xd} / (U_g - U_{xd})$

As $U_g \gg U_{xd}$: $R_x = R_h * U_{xd} / U_g$ (1)

Where $U_{xd} = (U_x - U_{ref})$ is voltage difference between **Ch** and **Ch_{ref}**.

Algorithm of the transhead impedance calculation

The next algorithm is implemented for impedance measurement and calculation:

1. Voltage ($U_x - U_{ref}$) is measured and filtered by 28-32Hz (-3dB) pass-band filter (8 order IIR Butterworth);
2. SQR value of ($U_x - U_{ref}$) is calculated for a 1 sec interval;
3. Impedance **R_x** is calculated once per second according formula (1).

This is repeated for all channels.

